medica+

JAPAN TECHNOLOGY

**Perfection Manual Massager**

**MassPro 5.0**

medica+

JAPAN TECHNOLOGY

# **Instruction Manual**


**Manual Massager  
MassPRO 5.0**

# Content

---


Packing list .....	01
User Instructions .....	02
Warnings .....	05
Declaration on Hazardous Substances ..	06
Warranty claims .....	07
Warranty terms .....	08

## Packing list

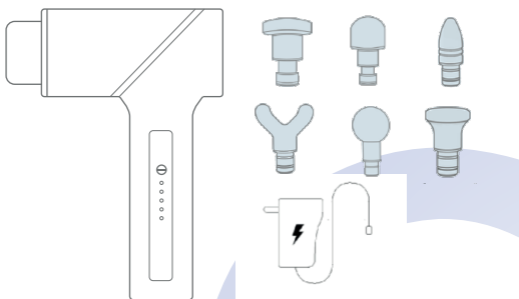
Model	MassPRO 5.0
Battery Capacity	3400 mAh
Power	15-60 W
Frequency	20 Hz, 26 Hz, 36 Hz, 46 Hz, 53 Hz
Battery Endurance	1.5-5.5 h
Torque Pressure	0-7.5 kg/cm
Amplitude	11 mm
Intelligent mode	10 min smart timing

# User Instructions
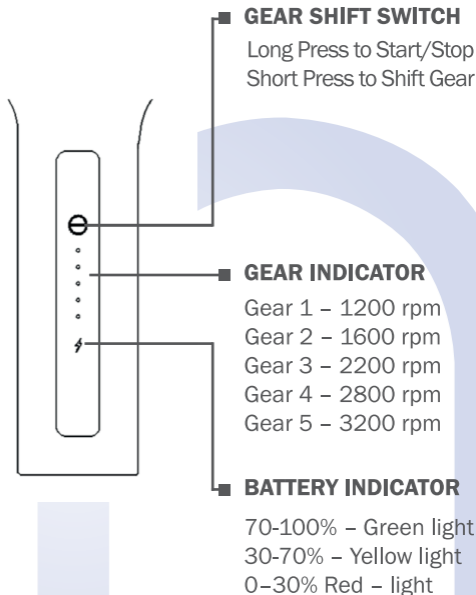
# User Instructions
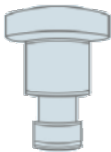
## **STAINLESS STEEL HEAD**

Strong massage and good feeling


## **BALL HEAD**

For large muscle groups (quads, glutes, etc) Or extream soreness

## **CONE HEAD**

Apply to strike the deep tissue, joint

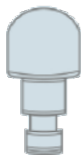
## **STAINLESS STEEL CYLINDER HEAD**

Strong massage  
and good feeling

## **FORK HEAD**

For neck, spinal,  
calves, etc
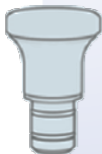
## **FLAT HEAD**

For general use, relax for  
whole body muscles

## Warnings

- Please use the original charger to charge the device, 3-5 hours is applied for a full charge, the indicator will be green when full charged.
- Do not use this machine while charging.
- Do not insert or touch the gap between the massage head and the shell with your fingers or other objects during use, so as not to pinch.
- This machine can not be self-disassembled or modified
- Due to the high power of the device, it is recommended to use it for 10~20 minutes at a time to avoid physical discomfort.
- DO not allow children to use alone.
- If you experience pain or discomfort be sure to immediately stop use.

### **CONSULT YOUR PHYSICIAN PRIOR TO USE IF THE FOLLOWING APPLY**

Osteoporosis patients, anyone with heart problems wearing an electronic medical device, pregnant women, anyone with physical trauma or any other medical complications must consult their physician prior to use.

# Declaration on Hazardous Substances

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## **THIS TABLE IS COMPILED ACCORDING TO SJ/T 11364.**

O: It means that the content of the restricted substance in all materials of part is less than the limit defined in GB/T 26572.

X: It means that the content of the restricted substance in at least on homogenous material of part is not less than the limit defined in GB/T26572.

# Warranty claims

## Warranty record

User name \_\_\_\_\_

Product number \_\_\_\_\_

Name model \_\_\_\_\_

Date of purchase \_\_\_\_\_

Date	Repair	Salesman	Remarks

## Warranty terms

1. If failure occurs due to manufacturer defect within 3 months of purchase date, the product must be complete and scratch free including all accessories and packaging. Products subject to inspection, if deemed tampered with or misused in any way this will automatically void the warranty.
2. A one year manufacturer's warranty is included with purchase and covers massage gun motor and battery failure due to manufacturer defect. Subject to inspection.
3. Warranty is only valid under normal personal use.
4. All damage caused by improper use is outside the scope of the warranty.
5. This card must be provided at time of warranty.
6. The manufacturer reserves the right to decline warranty services if deemed to be the result of improper use or negligence.

**THIS CARD MUST BE PROVIDED ALONG WITH  
RECEIPT OF PURCHASE**

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