ENG / UKR / RU

**Lifting
massager**

Skin Lifting 6.0

medica+
JAPAN TECHNOLOGY

medica+

JAPAN TECHNOLOGY

ENG

medica+
JAPAN TECHNOLOGY

**Instruction manual
lifting massager**


Skin Lifting 6.0

Content

| | |
|-----------------------------------|-------|
| Massager overview..... | 01 |
| Description, Basic operation..... | 02 |
| Mode selection and gear..... | 03 |
| Charge, Clean..... | 04 |
| Store, Discard..... | 05 |
| Method of use..... | 06-07 |
| Warning..... | 08 |
| Precautions..... | 09 |
| Common problem..... | 10-11 |
| Product parameters, Tips..... | 12 |
| Tips..... | 13 |

Massager Overview

Dear users: Thank you for purchasing our products! In order to facilitate your normal use of the lifting massager, please read the instruction manual carefully and keep it properly.

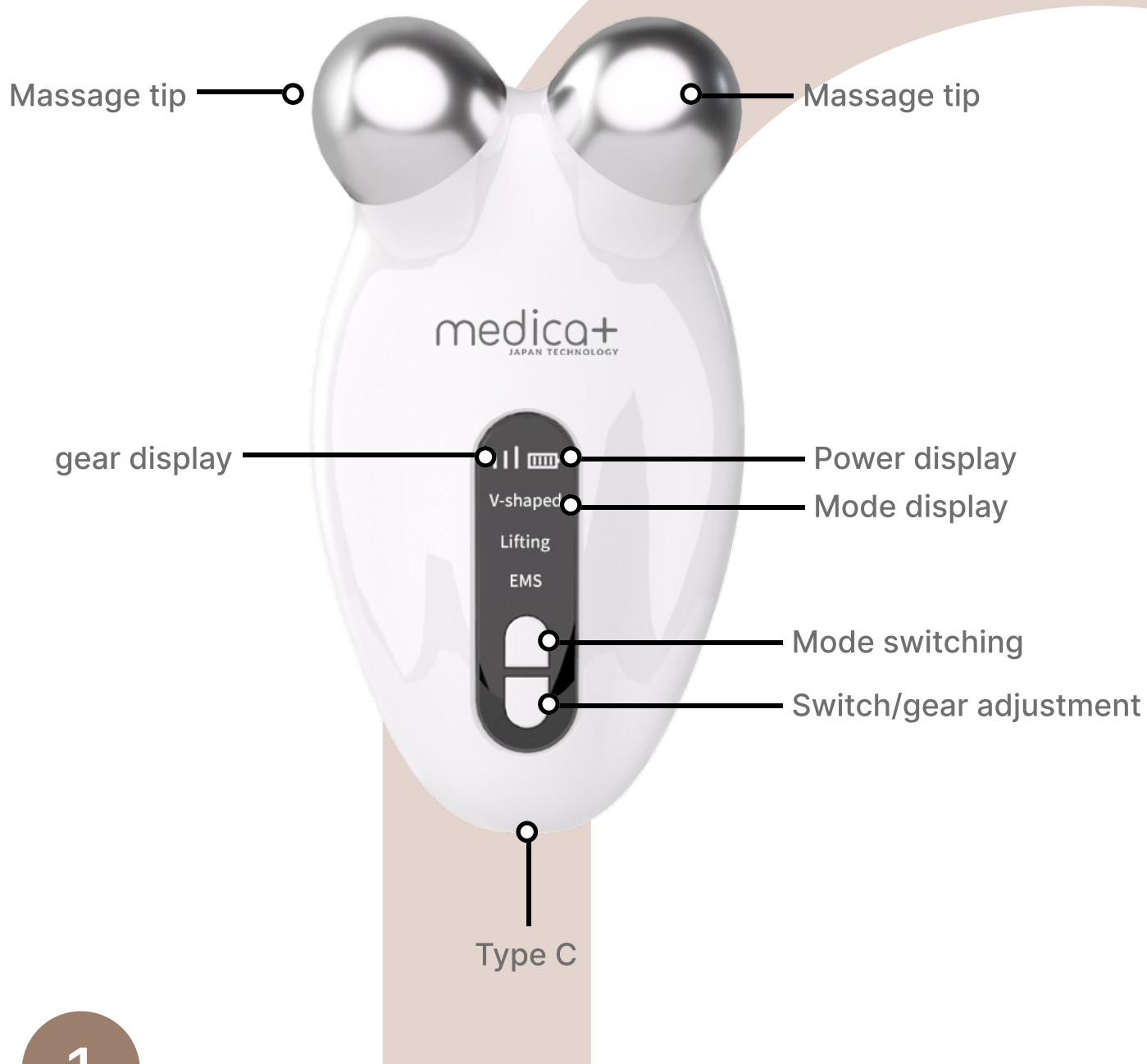Description

The microcurrent facial massager is a non-invasive and very safe home maintenance device. In just a few minutes, pulses of microcurrent are delivered deep into the skin to achieve unique results that can only be obtained in a beauty salon.

The special sphere design can conduct micro-current optimally. According to the skin feeling during use, you can also adjust the gear size at any time.

Basic operation

Power on/off method

- Press and hold the "⏻" button for 1.5 seconds to turn on the machine.
- Long press the "⏻" key for 1.5 seconds to shut down.

Mode selection and gear

1. The default "V-shape lifting" mode when starting up, micro-frequency vibration lifting and shaping;
2. Press the " M" button to switch to the "heat sense lifting" mode, constant temperature vibration lifting;
3. Press the "M" button again to switch to the "Lighten and Tighten" mode, and the EMS bionic microelectricity will tighten and lighten the lines.

The default is 1 gear when starting up, click the "U" button to switch 2 gears/3 gears

Charge

Use the charging cable distributed in the box to connect the charging hole at the bottom of the fuselage and the power adapter (need to bring your own) and turn on the power. The battery icon on the screen flashes to indicate charging. It takes about 1-2 hours to fully charge, and the battery icon is off. Blinking again means charging is complete

Clean

After use, just wipe the dirt on the body with a clean soft cloth, do not wash with water.

Store


When not in use, put it in a dry place out of the reach of children, and take it out and charge it every 1 month to maintain the normal performance of the battery.

Discard


This product contains a rechargeable battery. When the product is damaged or you decide to discard the product, according to your local regulations, you may need to take out the battery and send it to a designated place for unified disposal.

To ensure safety, please make sure that the battery is completely exhausted before removing it.


Method of use

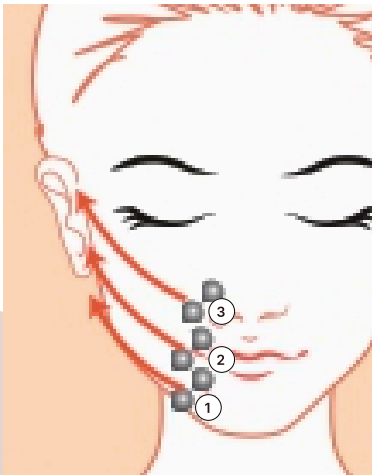Improve
nasolabial folds

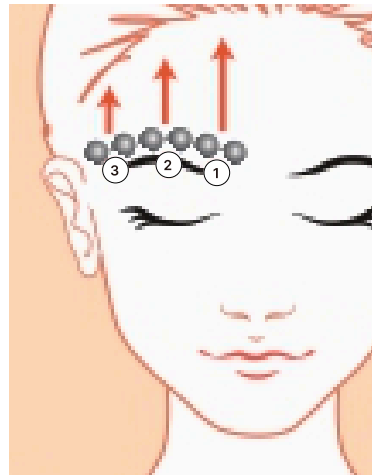Improve forehead
lines

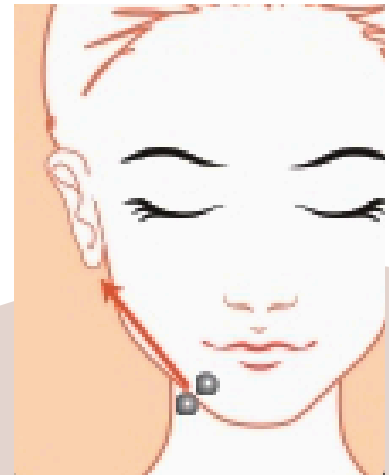Tighten the
contour

Improve nasolabial folds:

Place the ball of the massager on the cheeks ①, ②, and ③, and gently push in the direction of the arrow to the bottom, middle, and top of the ear, and repeat 3 times for 5 seconds each time.

Improve forehead lines:

Place the ball of the massager on the position ①, ②, and ③, above the brow bone, gently push along the direction of the arrow to the hairline, and repeat 3 times. Each time lasts 5 seconds.

Tighten facial contours:

Place the ball on the jaw line, pull it up and hold the ball for 5 to 10 seconds, repeat three times; the pressing and staying technique can also be used on other parts of the face that need firming and lifting.

WARNING: Do not use in:

1. At the midline of the neck (bones of the neck and throat area, possible effects on the thyroid gland)
2. Lip area and nasolabial folds (philtrum)
3. Edge around the eyes
4. Chest (used in the chest may affect the heart)
5. Groin area
6. Temple

Precautions

 The following people should not use this instrument:

1. Children
2. Pregnant women
3. Patients with epilepsy
4. Tumor and cancer patients
5. With a pacemaker
6. Implantation of a Cardiac Defibrillator
7. Others with electronic devices implanted in the body.

If you have any medical doubts, such as serious illness, epilepsy, etc., please consult your physician before use.

 The body, charger, and conductive massage head are strictly prohibited from entering water. It is recommended not to put it in the bathroom to avoid water ingress.

Common problem

Q: How long will it take to be effective?

A: Most customers can see the effect when they use it for the first time.

However, there are differences in individual conditions, and the effect will gradually appear after consistent use.

Q: How often is it recommended to do instrumentation?

A: In the first three months, it is recommended to use 5 times a week. After three months, you can use it 2-3 times a week to maintain the effect.

Q: How much current will I feel during use?

A: Microcurrent is a very low current value, which is basically not stimulating. Soothing and relaxing. A tiny tingling may be felt in some skin areas.

Of course, the sensitivity of each person's skin is different, and the current gear can be adjusted according to the skin feeling during use.

Q: Can I use it if I plan to have cosmetic surgery.

A: Yes. Microcurrent is used before surgery to improve skin firmness, reduce puffiness, and keep skin in the best condition. This device can be used in conjunction with other beauty and skin care procedures.

Product parameters

| | |
|-------------------|------------------|
| Products type: | Lifting massager |
| Products name: | Skin Lifting 6.0 |
| Input voltage: | 5V ==-1A |
| Battery capacity: | 400mAh |
| Rated power: | 3W |
| Charging time: | about 2 hours |

Tips

1. The warranty period of this product is one year.
2. Use it according to the precautions in the manual, and repair it free of charge when there is a non-human failure during the warranty period. (Please submit the product warranty card when repairing)

3. If there are the following situations, even within the warranty period, you must pay for maintenance:

(1) Failure and damage caused by misuse.

(2) Failures and damages caused by improper maintenance or modification of self-disassembly.

(3) Failures and damages caused by the use of an external power supply (voltage, frequency) due to irresistible factors such as fire, earthquake, electric shock, or abnormal voltage.

4. It is forbidden to use liquids with high acidity and alkalinity and corrosiveness.

5. This card will not be issued separately, please keep it properly.

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