
**ENG / UA / PL**


**Lifting  
massager**

**Skin Lifting 7.0**

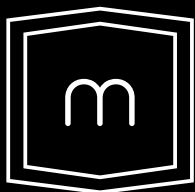**medica+**  
JAPAN TECHNOLOGY

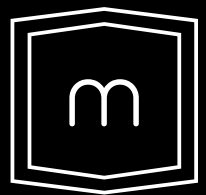medica+

JAPAN TECHNOLOGY

**ENG**

medicat+

JAPAN TECHNOLOGY

**Instruction manual  
lifting massager**

**Skin Lifting 7.0**

# **Contents**

---

Introduction .....	01-04
Names of Parts.....	05
Technical Specifications.....	06
Device Usage.....	07-12
Safety Instructions.....	13-14
Contraindications.....	15
Warranty Obligations .....	16
Package Contents .....	17

Dear Customer!

Thank you for choosing Medica+ products.

Please read this manual carefully before using the device.

## **Introduction**

This massager is specially designed by Medica+ for both home and professional use. It is intended for facial and neck massage. The device stimulates collagen and elastin production, provides lymphatic drainage and lifting effects, brightens pigmentation spots, regulates sebum production, minimizes pores, and improves skin tone and overall condition.

All stages of product development were carried out under the guidance of dermatology experts and engineers implementing advanced micro-vibration technologies. Thanks to its intelligent design and comprehensive approach, the Skin Lifting 7.0 microcurrent RF-Lifting facial massager provides highly effective skincare.

## **Principle of Acoustic Vibrations**

The Microcurrent RF-Lifting Massager generates up to 8,000 high-frequency vibrations per minute. This intense action helps deeply cleanse pores, removing impurities even from the deepest skin layers. As a result, the skin becomes smooth, fresh, and radiant, achieving a healthy appearance and natural glow.

## **Principle of Heat and Cold Alternation**

Repeated exposure to heat at 42°C stimulates collagen production, promoting skin cell regeneration, increasing elasticity, and helping to prevent wrinkles. The heat mode also aids in opening and cleansing pores.

Cooling down to 6°C instantly tightens pores, prevents moisture loss, and locks in nutrients within the skin. The alternation of warm and cold treatments enhances microcirculation in the deeper skin layers and accelerates cellular metabolism.

This skincare approach visibly reduces dark circles under the eyes and fine wrinkles around them.

The skin gains a natural glow and looks fresh even without makeup. Combined with micro-vibrations, the treatment enhances skin firmness and effectively relieves fatigue.

## **Principle of EMS (Electrical Muscle Stimulation)**


The EMS technology activates skin cells, stimulates collagen production, enhances elasticity, smooths wrinkles, and gives the skin a healthy, well-groomed appearance.

Microcurrents relieve facial muscle tension, reducing the number and depth of wrinkles while improving overall skin condition. At the same time, they stimulate the muscles that support the facial contour, providing a noticeable lifting effect in the treated area. Additionally, microcurrents boost elastin and collagen production, enhance lymphatic drainage and blood circulation, and accelerate cellular metabolism. This process improves nutrient exchange and speeds up skin regeneration, resulting in a healthier, more youthful appearance.

## Principle of LED Photon Therapy

- Blue light (470-480 nm) - penetrates to a depth of 0.25 mm, destroying harmful bacteria and preventing their re-growth. It has a strong anti-inflammatory and antibacterial effect. It helps to eliminate skin dullness and is effective in treating acne and post-acne.
- Red light (620-630 nm) - deeply penetrates the skin thanks to 18 LEDs, stimulating cellular metabolism, activating collagen production, improving blood circulation, and enhancing the effectiveness of cosmetic treatments. It tightens pores, regulates sebum production, increases skin elasticity, and is used in the treatment of acne.

# Names of Parts

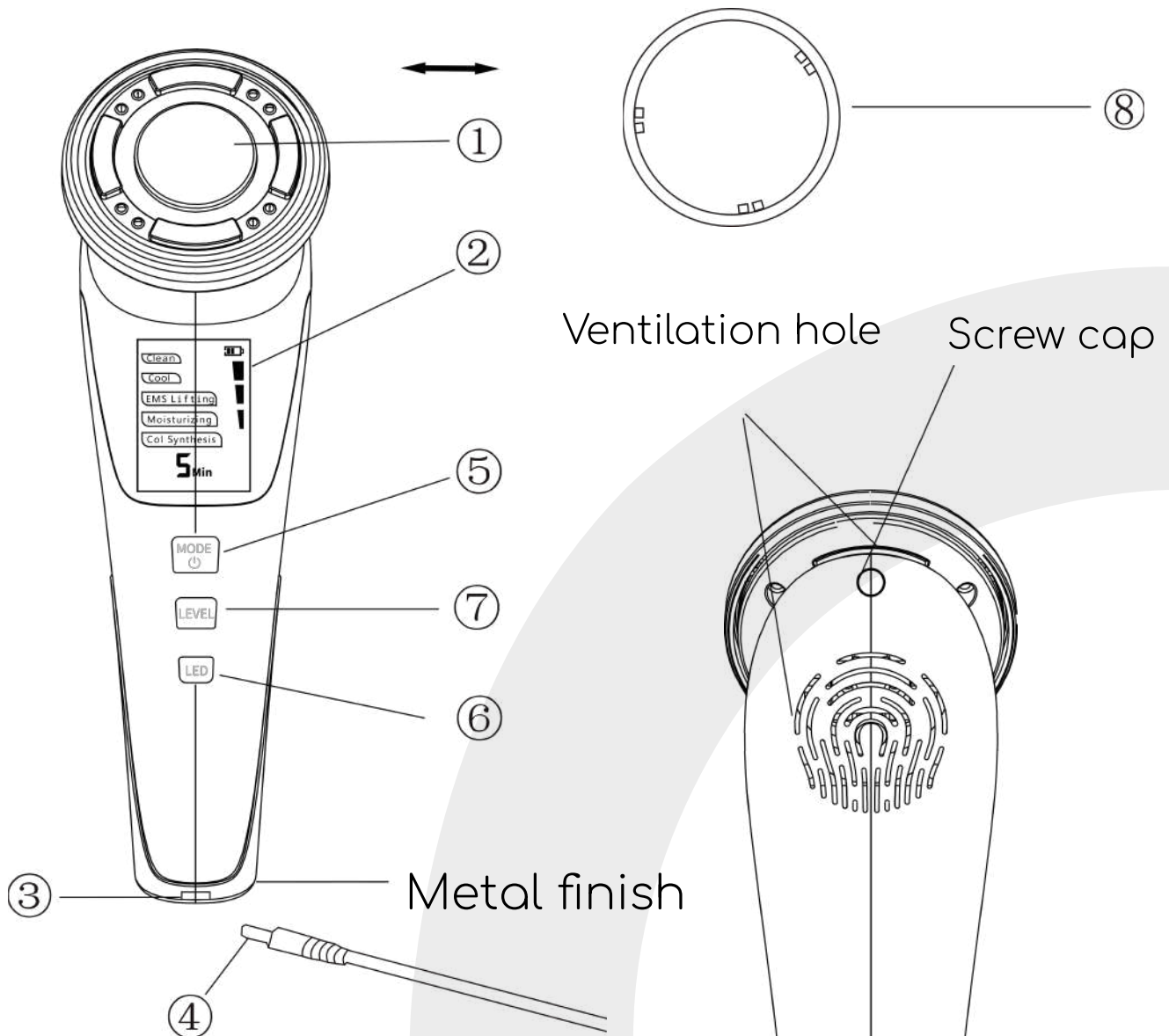1. Hot and cold panel
2. LCD display
3. Charging port
4. USB charging cable
5. Mode switch
6. LED therapy activation button
7. Power level adjustment
8. Sponge attachment ring

# Technical Specifications

**Name:** Skin Care Device

**Model:** Skin Lifting 7.0

**Power Supply:** 5V DC

**Temperature:** 6-42 °C ( $\pm 5$  °C)

**Vibration:** 8000 ( $\pm 10$ ) RPM

**Charging Time:** 2.5 hours

**Charging Current:** 500mA

**Power:** Up to 8.5W

**Battery Capacity:** 1000mAh

**Standby Mode:** 45 days

**Net Weight:** 145 g

**Usage Time:** 5 min/session or 2 min/session

**Device Size:** 175\*55\*67mm

**Box Size:** 215\*150\*80mm

# Device Usage


Before first use, fully charge the device. Charging time is approximately 2.5 hours. To turn on, press and hold the power button (5) for 3 seconds. The device will automatically start in the "Clean" mode at the first power level, using vibration and heating. The temperature can be adjusted according to the skin's needs by pressing the button (7). To change the LED light color, press the button (6). The "Clean" mode is equipped with an intelligent timer that automatically ends the procedure after 5 minutes.

## Cool Mode

Press the power button (5) again to activate the "Cool" mode. In this mode, the device uses vibration and cools the skin. The temperature can be adjusted using button (7). The intelligent timer automatically turns off the device after 2 minutes. After the procedure, the fan will continue to run for an additional 2 minutes to complete the cooling process.

## EMS Lifting Mode (Electromyostimulation)

Press the power button (5) to activate the "EMS Lifting" mode. The device generates EMS microcurrents and vibration. The intensity of the microcurrents can be adjusted using button (7). The mode is equipped with an intelligent timer for 5 minutes.

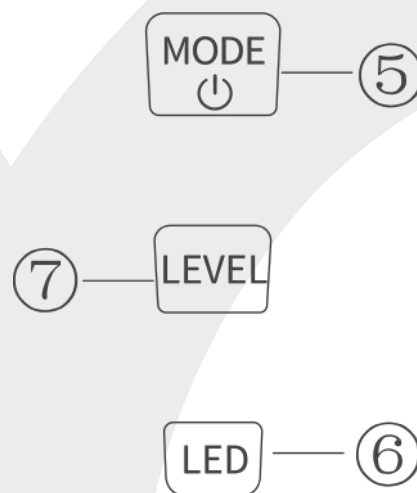
## Moisturizing Mode

Press the power button (5) to activate the "Moisturizing" mode. The device operates with vibration, heating, and EMS microcurrents. The intensity can be adjusted using button (7). The session duration is 5 minutes (intelligent timer).

## Collagen Synthesis Mode


Press the power button (5) to activate the "Collagen Synthesis" mode. In this mode, the device vibrates, cools the skin, and generates EMS microcurrents. The intensity can be adjusted using button (7). The duration is 2 minutes (intelligent timer).

## Instructions for Modes:

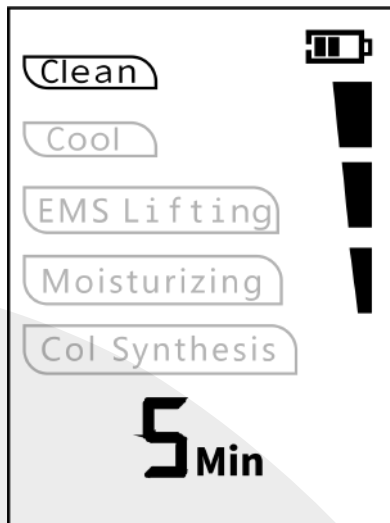
### 1st Mode "Clean" - Cleansing Mode

Used for additional skin cleansing. Moisten a cotton pad with toner. Secure the pad under the fixation ring. Perform the cleansing procedure, moving along the massage lines.

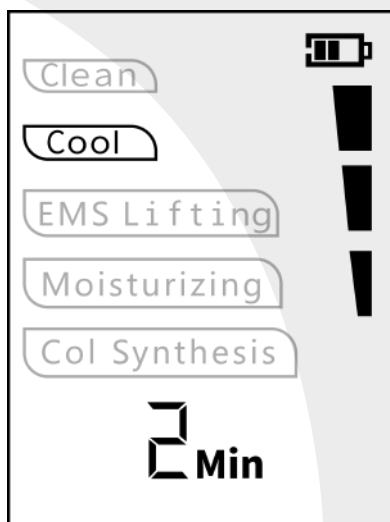
### 2nd Mode "Cool" - Cooling Mode

Cools the skin down to 6°C. Used for reducing puffiness, closing pores, preserving skin hydration, and toning the skin. It is used with Medica+ conductive gel, which is specially designed for this device to ensure maximum effectiveness.
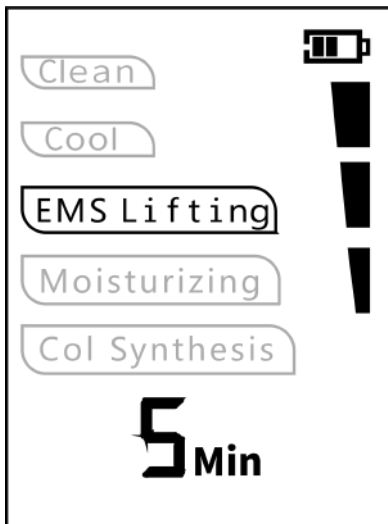
### 3rd Mode "EMS Lifting" - Electromyostimulation Mode

Performed with Medica+ conductive gel. This mode provides a pronounced lifting effect in the treated area, reduces signs of ptosis, improves lymphatic drainage and blood circulation, stimulates muscles, and smooths wrinkles.
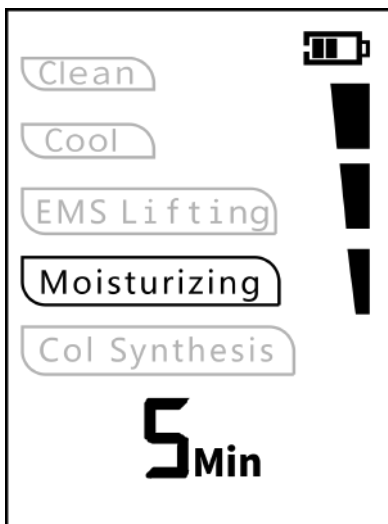
### 4th Mode "Moisturizing" - Moisturizing Mode

Heats the dermis to 42°C. Promotes deeper and more effective penetration of active ingredients from skincare products. Use after applying serum, fluid, or active ingredients suitable for your skin type.

Clean

Cool

EMS Lifting

Moisturizing


Col Synthesis

2 Min

## 5th Mode "Collagen Synthesis" - Collagen Synthesis Mode

Cools the skin down to 6°C. Helps to close pores, stimulates the production of collagen and elastin, and provides a rejuvenating effect.
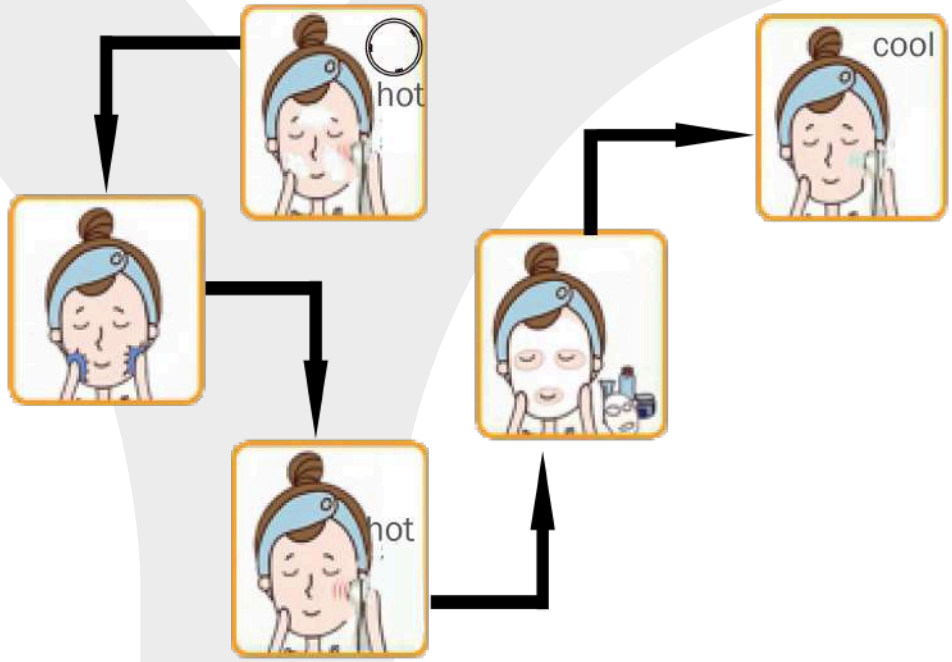
# Safety Instructions

- Please read this instruction before use. Follow the operating rules for safe use.
- Do not disassemble or modify the device yourself. In case of malfunction, contact the service center.
- Use the device only for its intended purpose according to the instruction.
- Avoid storing and using the device in high-temperature or high-humidity conditions (under direct sunlight, near an open flame, or in the bathroom).
- Do not use the device in wet conditions to prevent short circuits or electric shock.
- Keep the device out of reach of children and persons with limited capabilities. Children under 18 years old may use it only under adult supervision.
- After using the cooling mode, allow the device to cool down for at least 10 minutes before using it again.

## Warning:

- If the device overheats, immediately stop using it to avoid burns.
- Do not subject the device to impacts or throw it to prevent damage.
- Use recommended cosmetic products (serums, gels) for better results.
- Do not use the device if the power cord is damaged.

# Contraindications

- Pregnancy
- Cardiovascular disorders
- Arrhythmia
- Presence of a pacemaker, artificial heart, or lungs
- Presence of cancerous diseases
- Pathologies of internal organs
- Fever
- Acute phase of illness
- Increased sensitivity to skin diseases
- Sunburn
- Recently performed chemical peel, botox injections

# Warranty Obligations

Dear Customers!

Please keep this instruction carefully to use the warranty card. The warranty is provided for one year from the date of purchase (or for 2 or 3 years with the purchase of an extended warranty) in case of malfunction under normal operating conditions.

The warranty does not cover the following cases:

1. If the product is damaged due to improper use or lack of proper maintenance.
2. If the product is damaged as a result of disassembly or modifications.
3. If non-company parts were applied and used with the device.
4. Damage caused by force majeure circumstances.

According to the terms, technical services outside the warranty period are subject to additional charges.

Parts of the product replaced during the warranty period remain the property of MEDICA+.

# Package Contents

- Microcurrent RF Lifting Facial Massager Skin Lifting 7.0 – 1 unit
- Massager stand – 1 unit
- USB charging cable – 1 unit

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