**ENG / UA / PL**

**Manual  
Massager**


**MassHand Pro 6.0**

**medica+**

JAPAN TECHNOLOGY

**ENG**

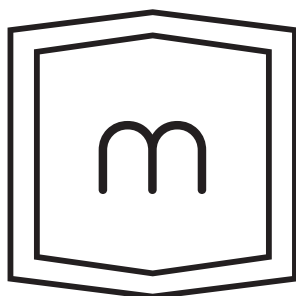**medicat+**  
JAPAN TECHNOLOGY

**Instruction**

**MassHand Pro 6.0**

# Content

---

<b>Precautions.....</b>	<b>01</b>
<b>Prohibition.....</b>	<b>01-02</b>
<b>Specifications.....</b>	<b>03</b>
<b>Product introduction.....</b>	<b>04</b>
<b>Massage heads .....</b>	<b>05</b>
<b>Instruction for use .....</b>	<b>06-07</b>
<b>Charging.....</b>	<b>07-08</b>
<b>FAQ.....</b>	<b>08</b>
<b>Troubleshooting.....</b>	<b>08-09</b>
<b>Maintenance .....</b>	<b>09</b>
<b>Enviromental protection.....</b>	<b>10-11</b>
<b>Warranty Terms .....</b>	<b>12</b>

Thank you for choosing our percussion massage gun. Before using it, please read this manual carefully to ensure proper operation and keep this manual for future reference.

- Relieving muscle soreness;
- Stimulating the circulation of blood;
- Relieving stiffness of cervical vertebra;
- Improving muscle recovery;

## **Precautions**

To avoid harm, fault, electric shock, fire, death and other equipment damage, please read the instructions and warnings carefully.

## **Prohibition**

1. Minors are forbidden to use this product. It is not a toy.
2. Forbid to put this product into water. If it falls into water accidentally, please send it for inspection in time.
3. Forbid to put this product near or into the fire.
4. Forbid to use it on head, neck, heart or genital.
5. Forbid to use it on bones. Use it on soft tissues only when not feeling pain or discomfort.
6. Forbid to use non-original charger to charge the product, otherwise it may cause a fire.

7. Forbid to drop the product intentionally.
8. Forbid to remove the screws or try to disassemble the product.
9. The product is not an alternative or complement of medical treatment. If you are feeling pain, undergoing medical treatment or using any other medical devices, please consult your doctor before using this product.
10. Forbid to put fingers, hair or other body parts near to the back of the massage heads in case of the squashing.

Do not touch this part when the machine is running without the massage heads being installed.



Do not touch this part when the machine is running with the massage heads installed.

# Specifications

Product size:	285*217*64mm
Color box size:	298*220*94mm (boxed)
Product weight:	800g/2000g (including packaging)
Speed:	2000-2800RPM/30 gears
Battery:	1500mAh/7.4V
Charging:	Type-C
Motor:	Brushed 12kg pressure bearing
Stroke:	6mm
Noise:	55dB
Body:	ABS/Metallic paint
Standard configuration:	Massage head 6

# Product introduction


Quick release massage head

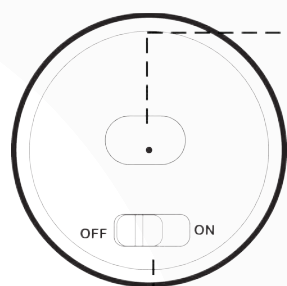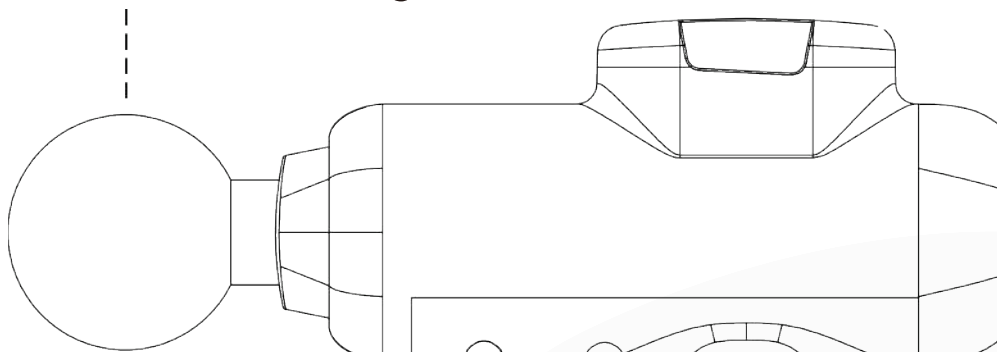
Charging port

Power switch


Silicone grip

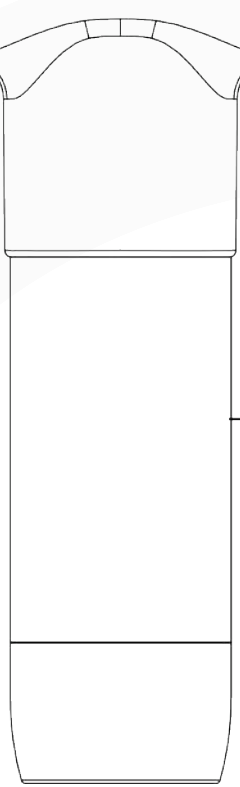Control panel

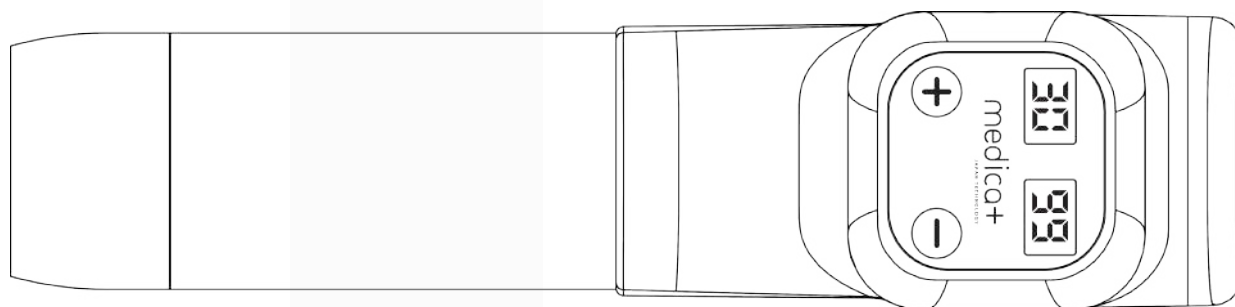
# Massage heads

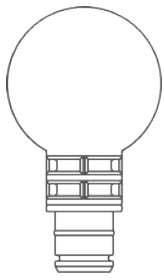
## Massage heads

For large muscle group massage and relaxation
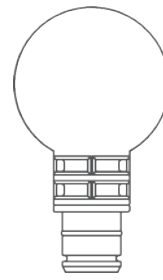
## Round head

For whole body muscle massage and relaxation
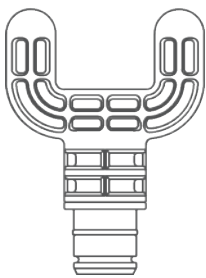
## Spinal head

For spine massage. Strong intensity

## Air cushion head

For soft tissue and whole body relaxation.

## Flat head


For large muscle group massage and relaxation

## Bullet head

For pain point relaxation. Strong intensity.

# Instruction for use
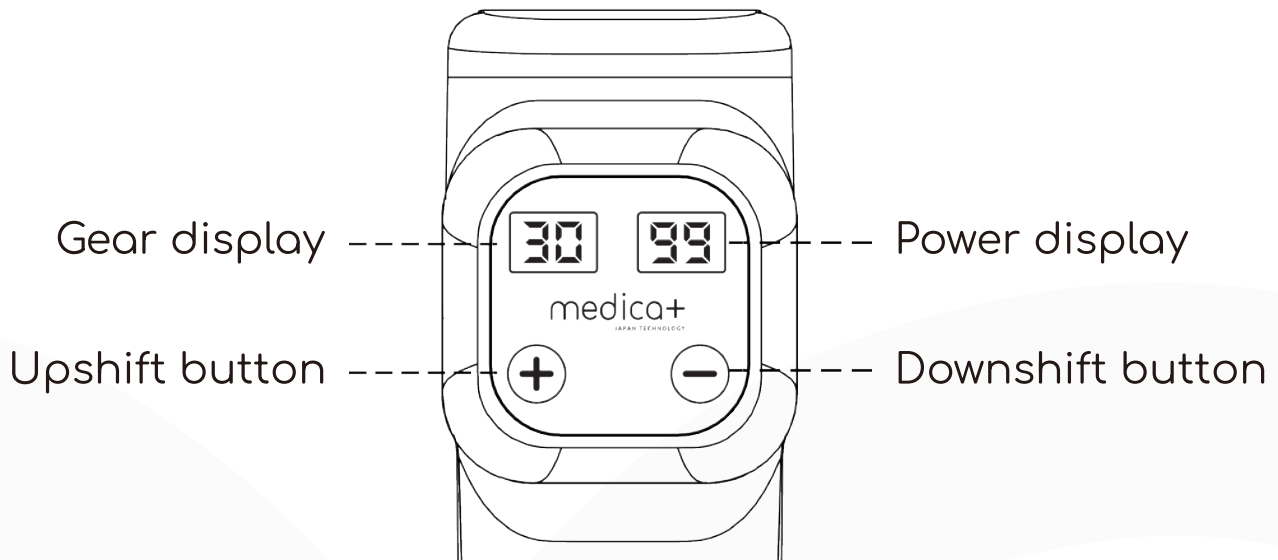
## Use of a massager

Use the "+" and "-" buttons to increase or decrease the vibration intensity with a short press on the gear button. Select the intensity that is comfortable for you and start the massage.

Getting the most out of your massager is easy if you follow a few simple steps:

1. Preparing for use: Select the desired massager attachment and insert it securely. This is the first step to an enjoyable massage.
2. Turning on the device: Press the "ON" button and your massager is ready to use. Two digits will appear on the display - on the right side the device's charge level, on the left side the vibration intensity mode (from 1 to 30).

3. Intensity setting: Easily adjust the intensity of the vibrations by pressing the gear buttons, and using the "+" or "-" keys to decrease or increase the intensity of the vibrations. Find the perfect strength that works for you for maximum comfort.
4. Start the massage: Start the massage by avoiding direct pressure on the spine and intervertebral discs. Allow the massager to relax your body, relieving you of tension.
5. End of the massage: After completing the session, press the shift button to reduce the intensity to the minimum, then turn off the massager with the "OFF" button.
6. Cleaning and care: If you used oil or massage cream, remember to wipe the attachments with disinfectant wipes to keep them clean and hygienic.

## **Charging**

The charging port is on the bottom of the massage gun handle.

- Please fully charge the battery for the first use
- Turn it off when charging
- Connect the power supply. Charging indicators will flash in turn and it shows that the machine is charging
- When the power display turns to 100%, it

means that the battery is full.

- Do not move the machine when it is charging.

## **FAQ**

### **Low temperature:**

Keeping the massage gun below 0°C for a long time will reduce the activity of battery and prevent normal functioning.

The product cannot work normally while being power on.

Solutions: Lay the product for two hours at temperature above 0°C, or raise the machine temperature fast with a hair-drier.

### **Low battery:**

When the battery is low, or user does not charge it for a long time, the battery will enter the low battery protection mode.

The product will not work normally at this time when you start it again.

Solutions: Turn it off with charger connected. Make sure it is fully charged, and it can run normally.

## **Troubleshooting**

Slight abnormal noise during use

Slight noise is a normal phenomenon in the loader after a period of use, please send it for inspection loader after a period of use, please send it for inspection

### **Battery does not charge**

Replace the battery please, for the aging or poor contact of battery.

## **Maintenance**

- When you do not use the product for a long time, turn it off and hold electricity while safe-keeping it. Store it in a safe, cool and dry place.
- Charge the product every half a year if you do not use it for a long time. Holding electricity prevents battery damage due to stains.
- Please choose neutral detergent and soft cloth for daily maintenance. Do not use engine oil, organic solvents or other chemicals for cleaning.
- Please turn off the power after using the product, or before cleaning the product.
- The product should be well fixed during transportation Avoid violent vibration, collision, rolling and falling.

# Environmental protection

Do not discard this product with general household waste at the end of its service life. Put this product in the recycle bin designated by the government or send this product back to the manufacturer's after-sales service center. The service center staff will recycle it in an environmental-friendly manner.

## Name of harmful substances in the product

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This form is based on the rules of ST 11364

○: The content of harmful substances in all parts meet the limited requirements of GB/T 26572

×: The content of harmful substances in at least one part exceeds the limitation of GB/T 26572

## Environmental protection use period

The time limit of environmental pollution, serious damage to people and property, or the leakage of toxic and harmful substances, is 10 years in normal use as described in this manual.

This sign indicates that the environmental protection period of this product is 10 years in normal use. Some components may also have environmental protection period signs. The specific period of environmental protection for each component shall be subject to its own sign.

# Warranty Terms

Dear Customers!

Please keep this instruction properly to use the warranty certificate. The warranty is valid for one year from the date of purchase (or for 2 or 3 years if an extended warranty is purchased) in case of malfunction under normal operating conditions.


The warranty does not cover the following cases:

1. If the product is damaged due to improper use or lack of proper maintenance.
2. If the product is damaged due to disassembly or modifications.
3. If non-original parts that do not belong to our company have been applied or used in the device.
4. Damage caused by force majeure circumstances.

According to the terms, any technical maintenance provided outside the warranty period is subject to an additional charge.

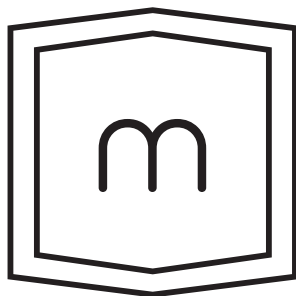Parts replaced during the warranty period remain the property of MEDICA+.

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