
ENG / UA / PL


Oral irrigator

ProWater Station 7.0

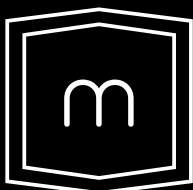
medica+
JAPAN TECHNOLOGY

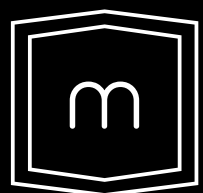medica+

JAPAN TECHNOLOGY

ENG

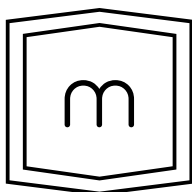medica+
JAPAN TECHNOLOGY

Instruction manual

ProWater Station 7.0

Contents

Package Contents	01
Safety Precautions	02-04
Schematic Diagrams	05-07
Usage Instructions	08-12
Recommended Technique	13-17
Cleaning and Maintenance	18
Warranty Terms	19
Common Malfunctions and Solutions	20
Technical Specifications	21

Dear Customer!

Thank you for choosing Medica+ products.

Please read this instruction manual carefully before using the device.

Package Contents:

Main unit with reservoir – 1 pc.

Instruction manual – 1 pc.

Classic jet tip – 3 pcs.

Tongue cleaner tip – 1 pc.

Orthodontic tip – 1 pc.

Plaque removal tip – 1 pc.

Periodontal gum tip – 1 pc.

Safety Precautions


When using the device, especially around children, it is important to follow basic safety measures, including:

- ✓ Always turn off the device immediately after use.
- ✓ Do not use it while bathing.
- ✓ Do not leave or store the device in places where it may fall into a bathtub or sink.
- ✓ Do not immerse the device in water or any other liquid.
- ✓ If the device falls into water, do not touch it while it is turned on – first, disconnect the power supply to reduce the risk of burns, electric shock, fire, or injury.
- ✓ Do not leave the device turned on unattended.
- ✓ Children, individuals with disabilities, or those with reduced self-care abilities should use the device only under supervision.
- ✓ Use the device only for its intended purpose, as described in the manual.


- ✓ Do not use the device if its cord or plug is damaged, if it malfunctions, has been dropped, damaged, or submerged in water. In such cases, contact a service center.
- ✓ Keep the cord away from heated surfaces.
- ✓ Do not cover ventilation openings or place the device on soft surfaces (such as a bed or sofa) to prevent airflow blockage. Do not allow hair or foreign objects to enter the openings.
- ✓ Do not use the device while sleeping or in a drowsy state.
- ✓ Do not insert or drop foreign objects into the device.
- ✓ Do not use the device outdoors or in areas where aerosol sprays are used or oxygen is being supplied.
- ✓ Connect the device only to a properly grounded outlet.
- ✓ Always turn off the device before filling the reservoir. Fill it only with water (a small amount of mouthwash may be added). Follow the instructions for proper filling levels.

- ✓ During initial use, minor gum bleeding may occur – this is a normal reaction depending on gum sensitivity.
- ✓ Children should not play with the device. It must be stored in a place inaccessible to children.
- ✓ If the power cord is damaged, it must be replaced at a service center.

Schematic Diagrams

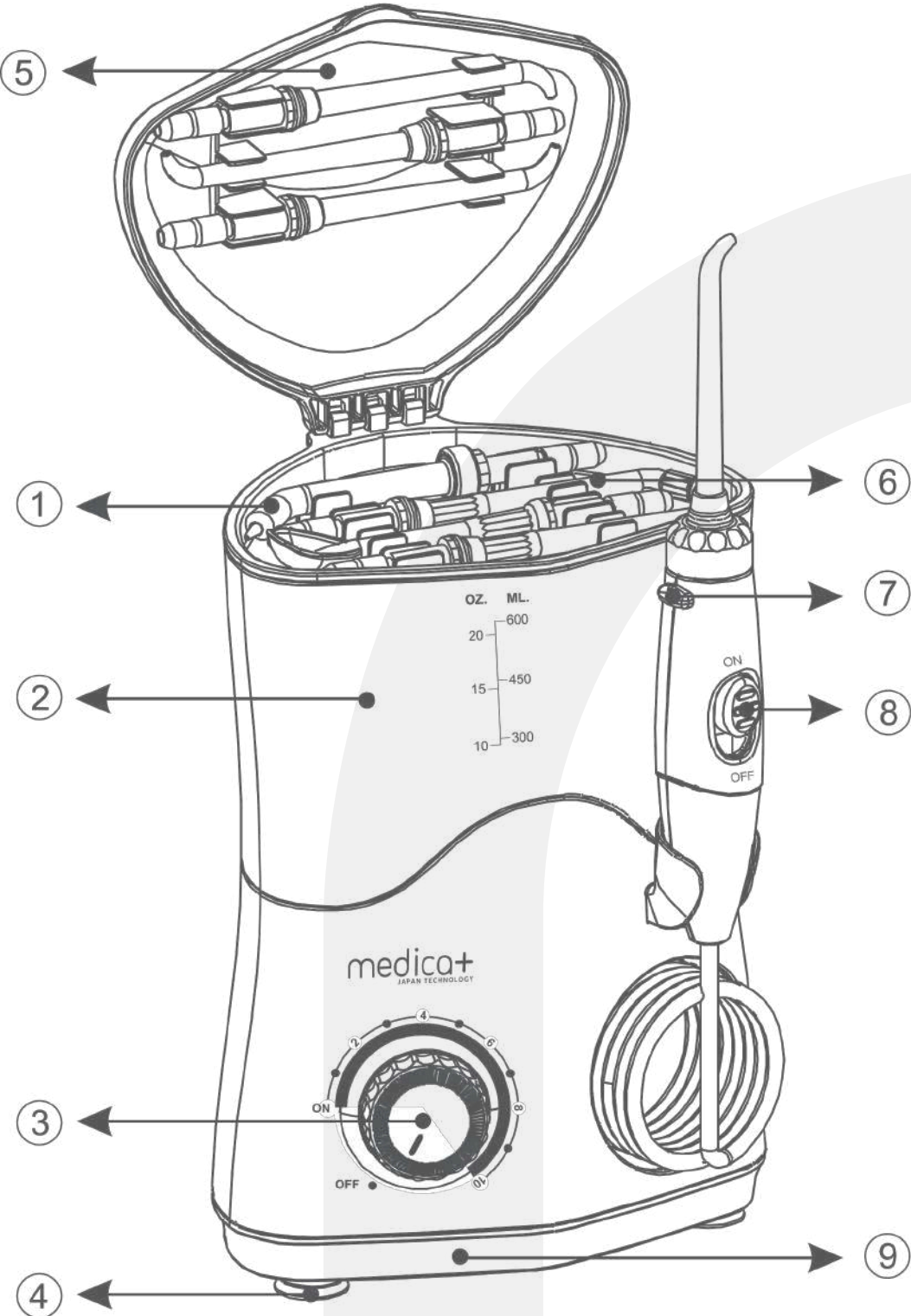
1. 7 nozzles
2. Reservoir
3. On/Off
4. Suction cup
5. Reservoir lid
6. Nozzle
7. Nozzle holder handle
8. Pause
9. Base suction cups

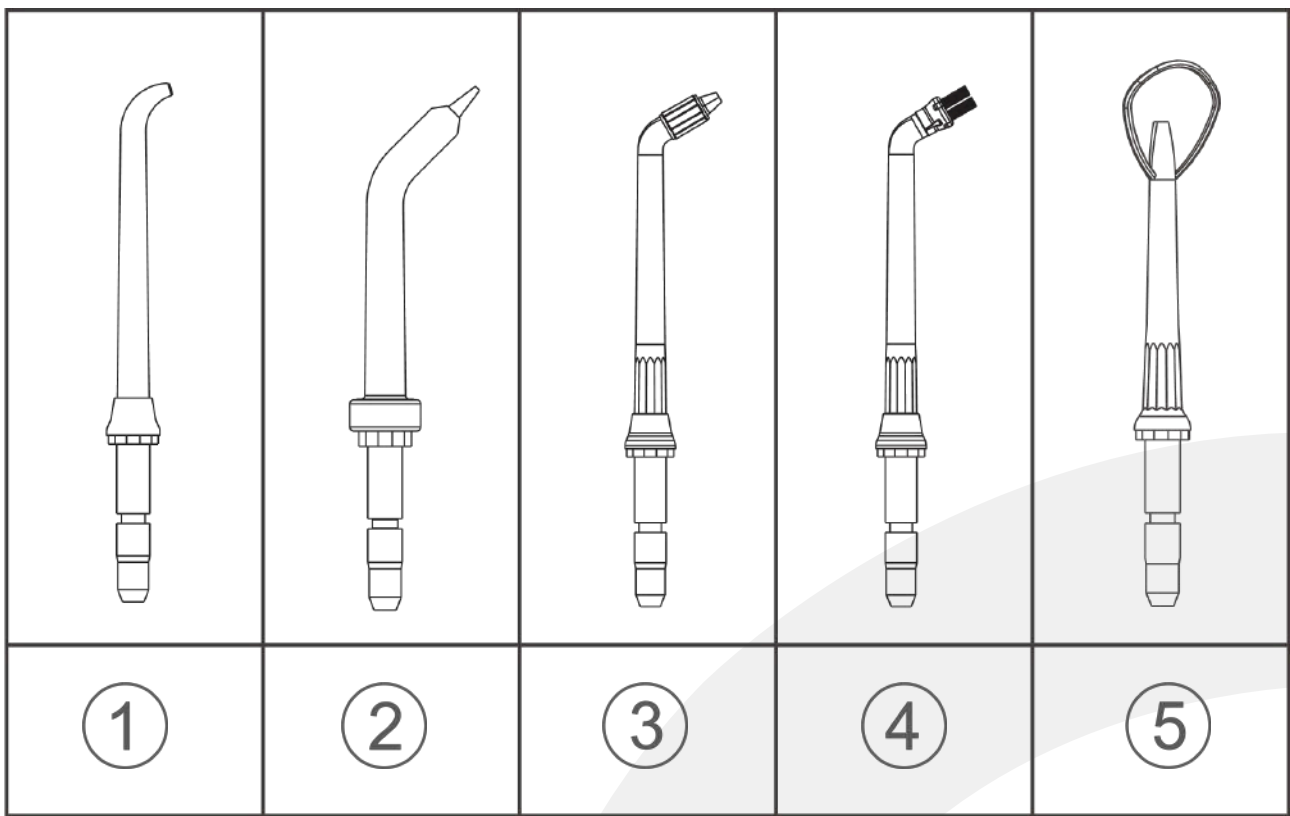1. Classic jet nozzle
2. Periodontal gum nozzle
3. Orthodontic nozzle
4. Plaque removal nozzle
5. Tongue cleaning nozzle


Usage Instructions

Plug the device into the power outlet. If the outlet is connected to a wall switch, make sure it is turned on.

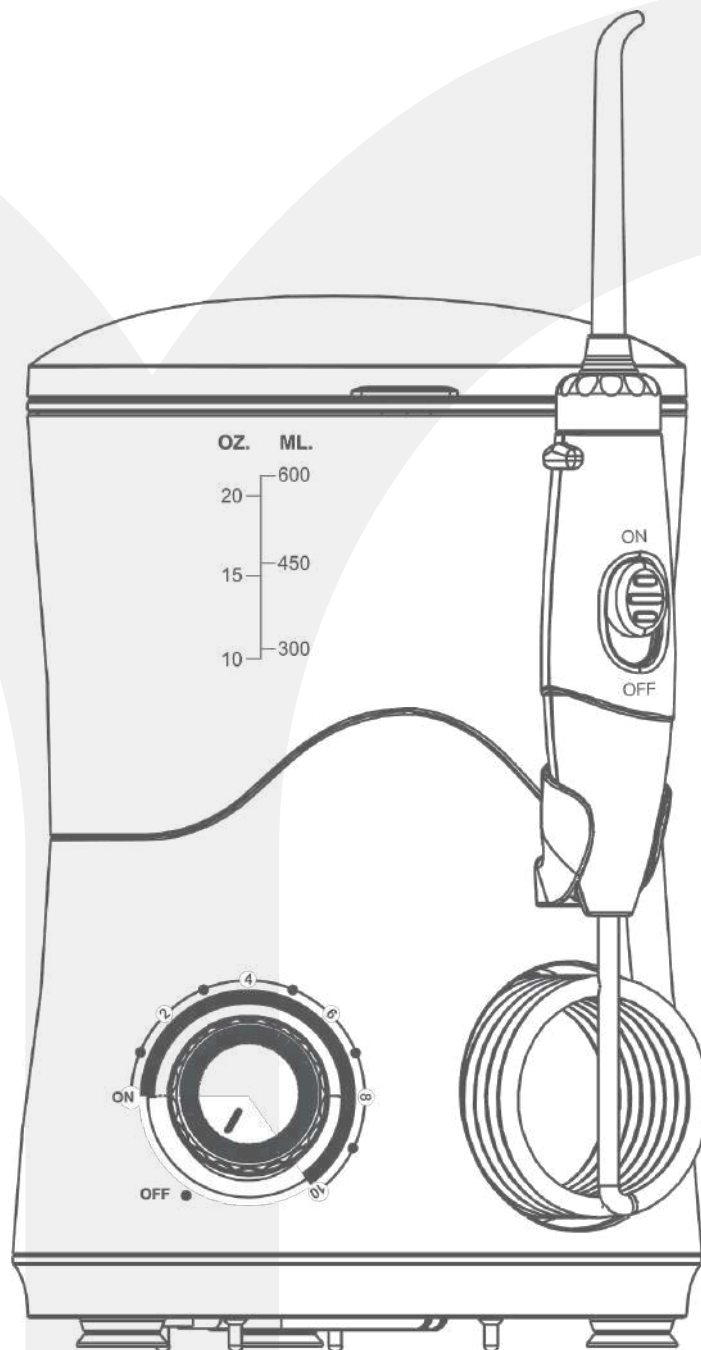
Preparation for Filling

Remove the lid and the reservoir from the base, fill the reservoir with warm water, place it back onto the base, press firmly, and put the lid on.

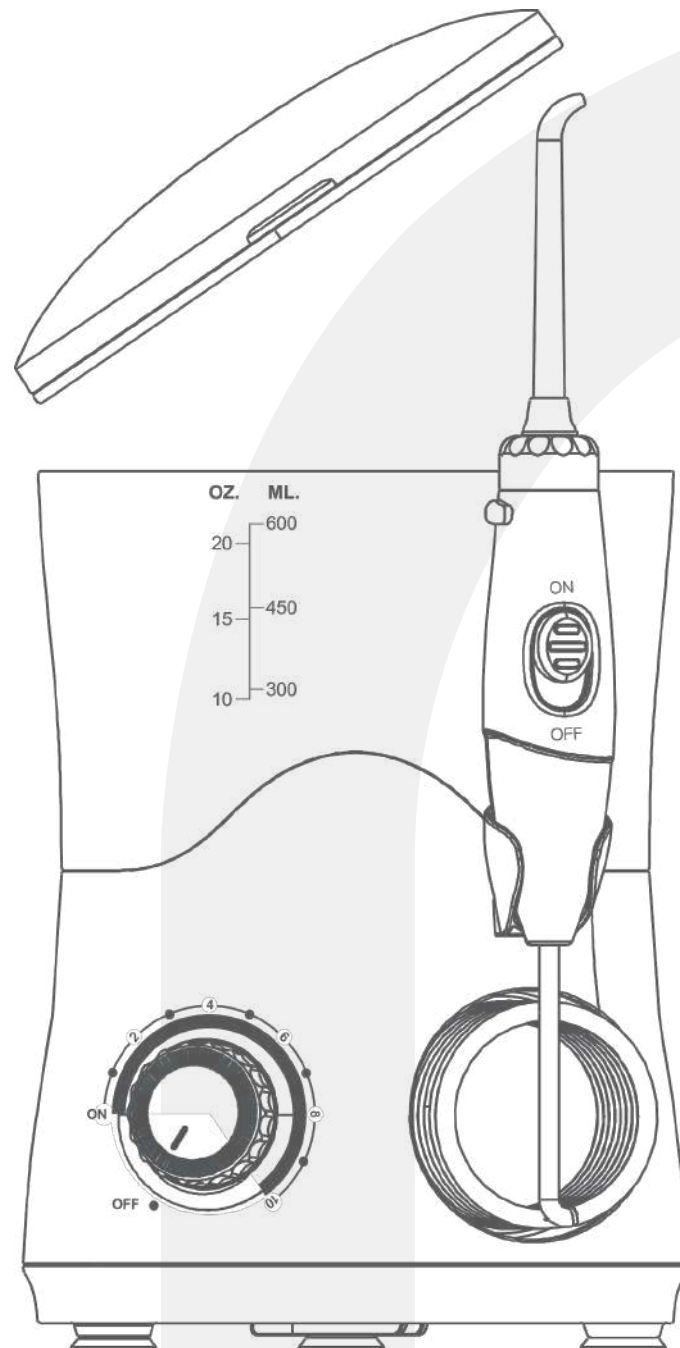
Installing and Removing Tips

Take a tip and insert it into the top of the device handle. Press until the tip clicks into place—the ring should align with the end of the handle. To remove, press the tip release button on the handle and gently pull the tip upward.

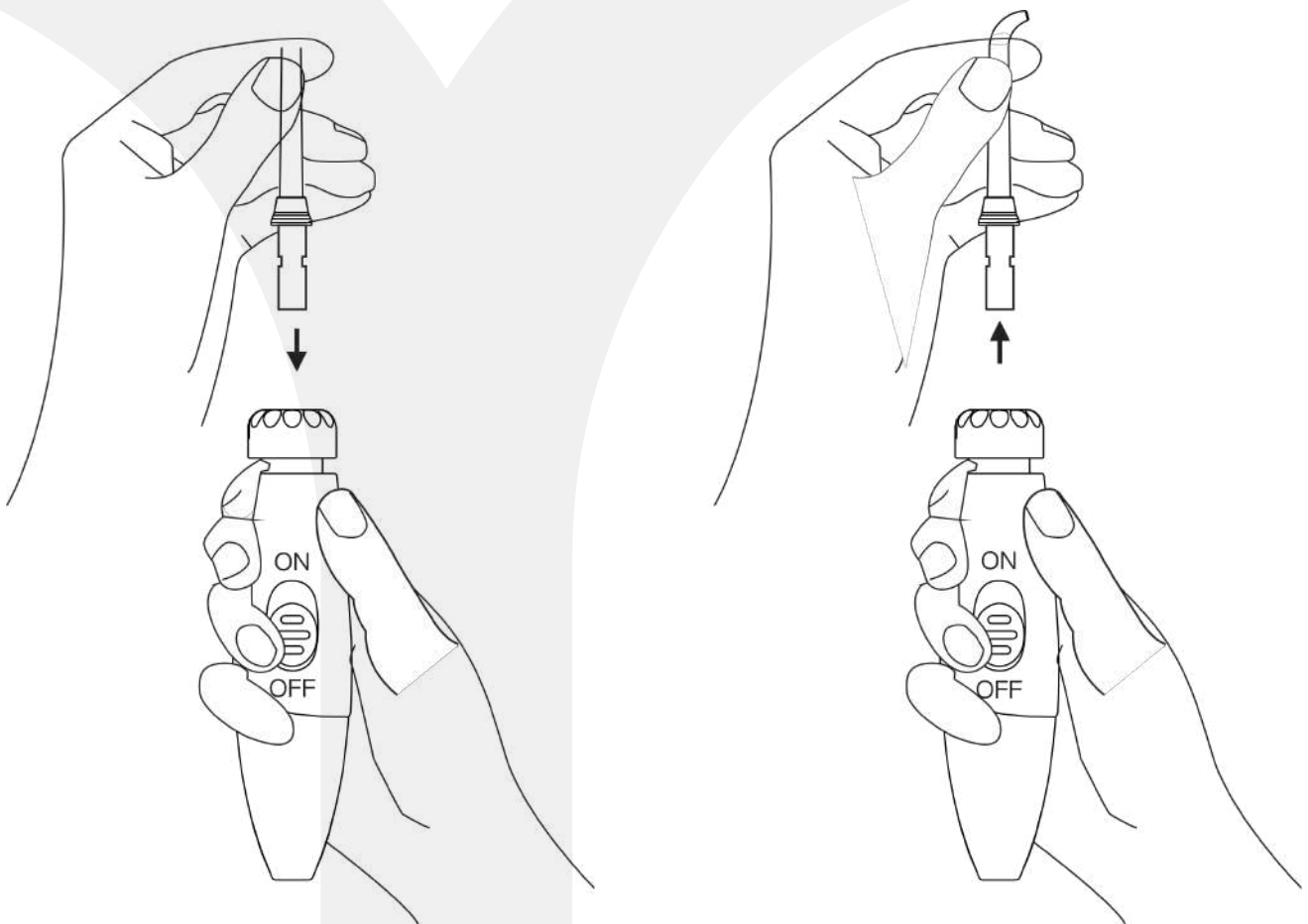
Pressure Adjustment

When you hear a "click," it means the switch is off. To adjust the pressure, turn the regulator clockwise—from the minimum to the maximum level. The device will automatically turn off after 3 minutes of inactivity. To turn it on again, rotate the regulator back to the "click" position.

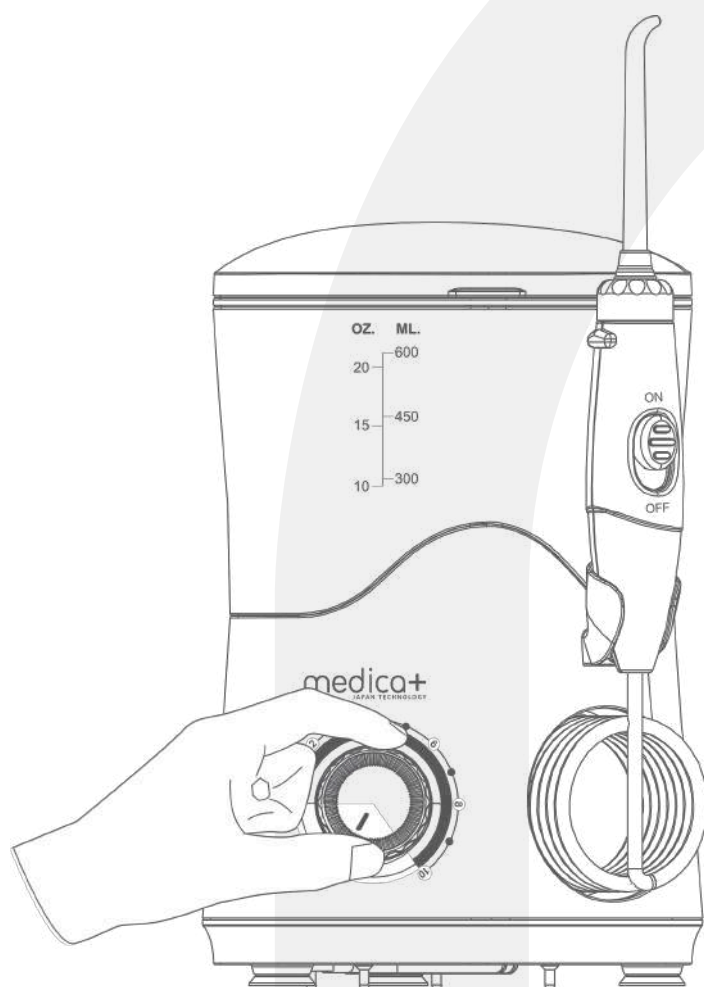
Using the Oral Irrigator

Lean over the sink, insert the nozzle into your mouth, and aim it at your teeth. Use your free hand to turn on the device and set the desired pressure level. To control the water flow, rotate the device handle.

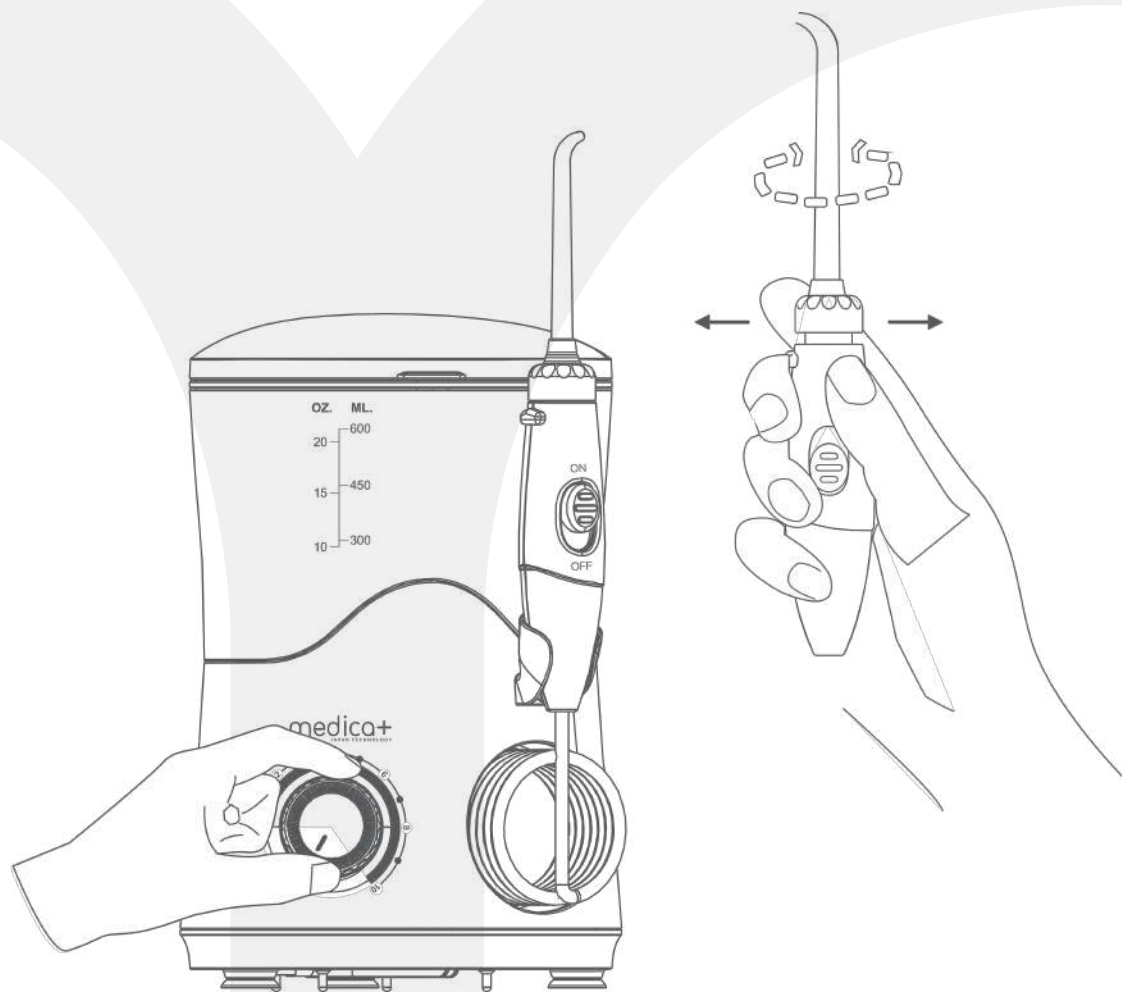
Recommended Technique

Direct the water jet at a 90° angle to the gum line, keeping your lips slightly closed to prevent splashing while allowing water to flow freely into the sink. Start cleaning from the back teeth, gradually moving toward the front. Slowly guide the nozzle along the gum line, making brief stops between teeth. Continue the procedure until all areas between and around the teeth are thoroughly cleaned.

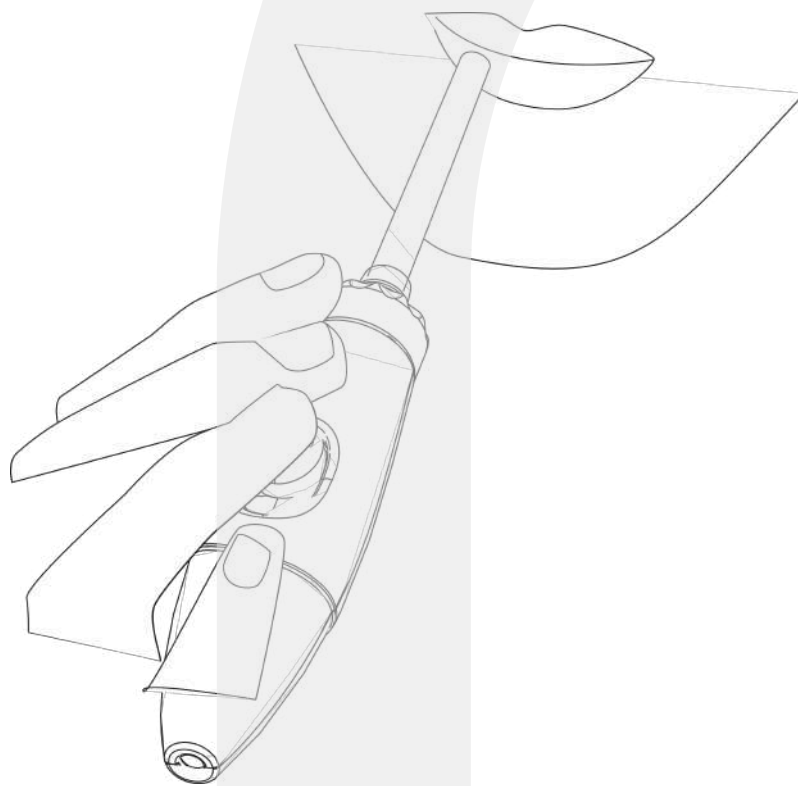Pause

You can temporarily stop the water flow at any time by sliding the button on the device handle downward.

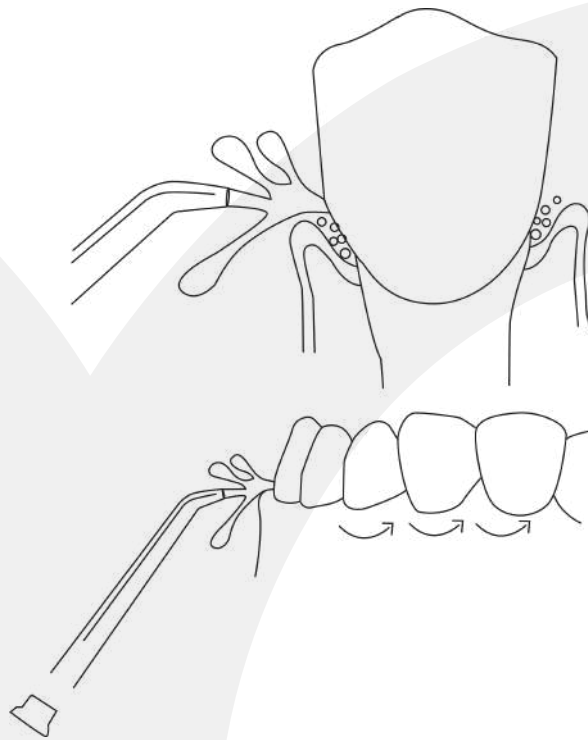Completion of the Procedure

Turn off the device. Drain the water from the reservoir and unplug the device from the socket.


Periodontal Tip

Allows directing water or liquid oral medications directly under the gumline, ensuring deep cleaning of hard-to-reach areas.


Classic Jet Tip

Cleans the oral cavity at a 360° angle, reaching the deepest areas, including tooth crevices and spaces under the gums. Removes food debris, bacteria, and plaque, providing a feeling of cleanliness.

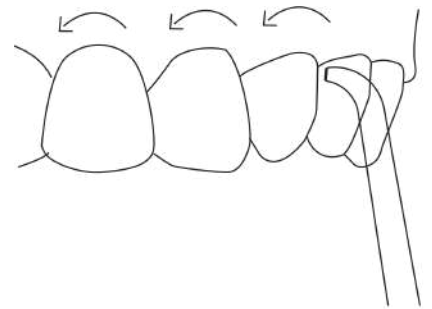
Tongue Cleaning Tip

Before turning on the irrigator, place the tongue scraper on the central back part of the tongue. Apply slight pressure and slowly move it forward to remove plaque and refresh breath.

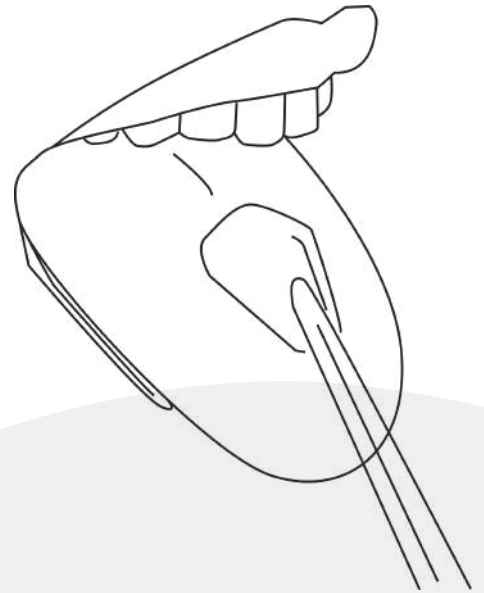
Plaque Removal Tip

Effectively combats soft plaque before it hardens. This specially designed tip cleans hard-to-reach areas, especially around braces and other orthodontic constructions.

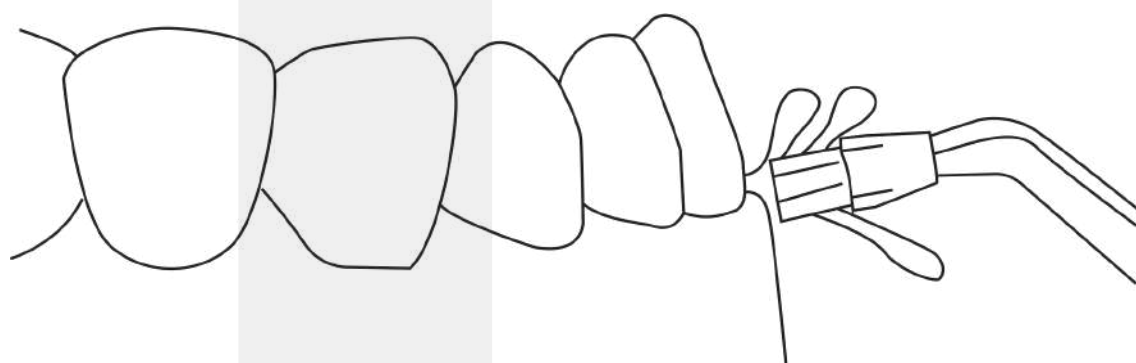Orthodontic Tip

Designed for cleaning teeth with braces and other dental appliances. Use the bristle brush to clean around the braces while allowing the water flow to wash away all debris.

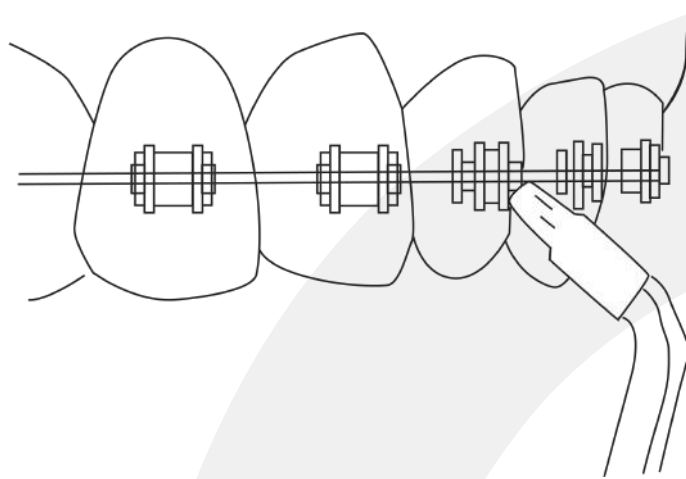Cleaning and Maintenance

After using special solutions, the irrigator must be rinsed. To do this, fill the reservoir with warm water, direct the stream into the sink, and turn on the device until all the water has been expelled. This will help prevent clogging and extend the device's lifespan. The housing can be wiped with a soft cloth and a non-abrasive cleaning agent. Always unplug the device before cleaning.

Warranty Terms

Dear Customers!

Please keep this manual properly to use your warranty card. The warranty is provided for one year from the date of purchase (or for 2 or 3 years if an extended warranty is purchased) in case of malfunction under normal operating conditions.

The warranty does not cover the following cases:

1. If the product is damaged due to improper use or lack of proper maintenance.
2. If the device has been damaged due to disassembly or modifications.
3. If non-original parts that do not belong to our company have been used.
4. Damage caused by force majeure circumstances.

According to the terms, technical maintenance provided outside the warranty period is subject to an additional fee.

Any parts replaced during the warranty period remain the property of MEDICA+.

Common Malfunctions and Solutions

Problem	Possible Cause	Solution
Leakage from the tip and tongue cleaner between the upper part and the handle.	The nozzle is not securely installed.	Remove and reinstall the nozzle.
The periodontal nozzle has become too soft.	The nozzle is worn out.	Replace it with a new one.
The irrigator is leaking	Hose damage	Can be replaced at the service center.
Leakage in the reservoir.	Missing water reservoir valve or the water reservoir valve is inverted.	Replace the reservoir valve at the service center.
Inadequate pressure	The reservoir is incorrectly positioned	Clean the filtration mesh
No water flow	The motor is not working	Make sure the outlet is working by testing it with another device (check the wall switch).

Technical Specifications

Weight – 0.930 kg

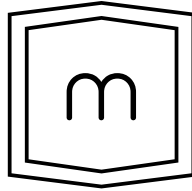
Pressure Range – 30-125 psi (pounds per square inch)

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