**ENG / UA / PL**

**Ultrasonic face  
scrubber**

**VibroSkin 8.0**

**medica+**  
JAPAN TECHNOLOGY

medica+  
JAPAN TECHNOLOGY

**ENG**

medica+  
JAPAN TECHNOLOGY

**User manual**

**VibroSkin 8.0**

# Content

---

Description.....	01-02
Important Safety Measures.....	03
Indications .....	04
Contraindications .....	05
Usage Recommendations .....	06
Mode Description .....	07-10
Technical Specifications.....	11
Warranty Terms .....	12
Package Contents.....	12

Dear Customer!

Thank you for choosing Medica+ products.

Please read this manual carefully before using the scrubber.

## **Description**

Ultrasonic vibrations are high-frequency mechanical waves with intense oscillations widely used in medicine and cosmetology. By generating a small amount of energy and high-frequency vibrations, these waves create a gentle resonant effect on the skin, improving cell activity and the permeability of fatty tissues while cleansing sweat glands and tightening pores.

The ultrasonic facial scrubber utilizes nanotechnology to generate vibrational oscillations at a frequency of 23,000 Hz per second, causing intense friction between molecules and accelerating microcirculation and metabolism in skin cells. The deep high-frequency vibration massage enhances microcirculation, promotes the breakdown of fat cells, and facilitates their removal through metabolic processes.

The spatula-shaped attachment ensures fast and effective pore cleansing, removing dirt, acne, dead skin cells, sebum deposits, and cosmetic residues.

Metal blade

On/Off button

Peeling Ion+ (Cleansing mode)

Peeling Ion- (Moisturizing mode)

Lifting Charge (Lifting mode)

Battery indicator

Protective cap

# Important Safety Measures

1. The device is not fully charged at the factory. Charge it before first use.
2. Before charging the device, ensure that the USB charging cable is not damaged or deformed. Do not connect the device to the USB charging cable with wet hands.
3. Before storing the device, make sure it is turned off. Store it in a dry place out of reach of children.
4. Use antiseptic wipes or a cotton pad moistened with a chlorhexidine solution to clean the body and spatula of the ultrasonic scrubber. Do not use corrosive cleaning agents. Clean the device after each use.
5. Do not throw the device into fire or store it in high-temperature or humid environments.
6. Do not immerse the device in water or rinse it directly under a faucet.
7. If the device has been damaged due to a fall, do not use it.
8. Do not disassemble or repair the product yourself. For repairs, please contact the Medica+ service center.

# Indications for Using the Ultrasonic Scrubber

- Dull facial skin tone
- Enlarged pores
- Oily skin
- Acne and comedones
- Scars and post-acne marks
- Fine wrinkles
- Loss of firmness
- Increased pigmentation

# Contraindications for Using the Ultrasonic Scrubber

Although the use of an ultrasonic scrubber is a relatively safe procedure, even for home use, it does have several contraindications. The device is not recommended for individuals with:

- Couperose (vascular skin conditions)
- Active herpes outbreaks
- Fever
- Malignant tumors
- Cardiovascular diseases
- Arterial hypertension
- Epilepsy
- Trigeminal neuralgia
- Pregnancy (at any stage)
- Increased skin sensitivity
- Sunburn or recent chemical peels
- Tendency to allergic skin reactions
- Neurological and dermatological conditions

# Usage Recommendations

Remove the protective cap before using the scrubber.

Press and hold the power button for 1.5 seconds to turn the device on/off.




Press the button briefly (1 second) to switch between modes and select the desired one (Peeling Ion+ / Peeling Ion- / Lifting Charge).

# Mode Description

## Mode 1 - Peeling Ion+ (Ultrasonic Cleansing)

This mode is used for deep facial cleansing – removal of blackheads, comedones, pore cleansing, exfoliation of dead skin cells, removal of excess oil and makeup residues, reduction of pigmentation spots and freckles.

In this mode, hold the scrubber with the sharp edge of the spatula against the skin at a 45-degree angle. Movement of the spatula – from the periphery to the center.

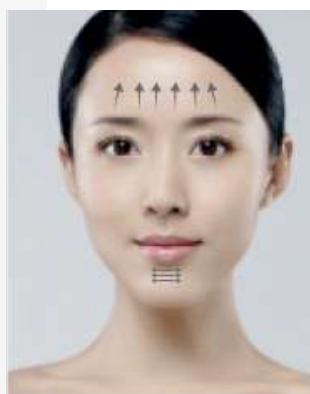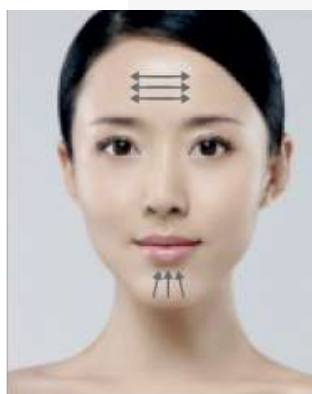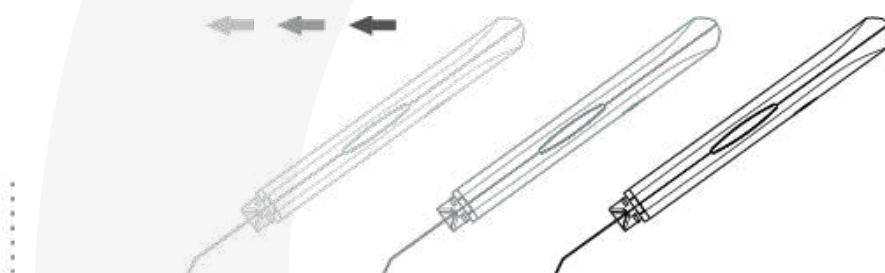Before the cleansing procedure, it is recommended to apply Medica+ Hydration Gel to the skin for 15 minutes under a film. After this period, remove the film locally and begin the cleansing procedure. If necessary, additionally moisturize the skin with a cotton pad soaked in toner or warm water.

After the procedure, apply Medica+ Soothing Mask or Medica+ Pore Control Mask to the skin.


- For oily and combination skin, the procedure can be performed 1-2 times a week.
- For dry skin, the procedure should be performed once a week.

Procedure duration is individual and depends on the degree of skin contamination. On average, it should not exceed 10-15 minutes.

## **Mode 2 - Peeling Ion- - Hydration Mode**

This mode is used for skin massage. High-frequency vibration improves cell activity, enhances blood circulation, and accelerates metabolism. As a result, facial skin becomes more elastic, smooth, and soft.

In this mode, hold the scrubber with the flat side of the spatula against the skin. Movement of the scrubber – slowly from the center to the periphery.

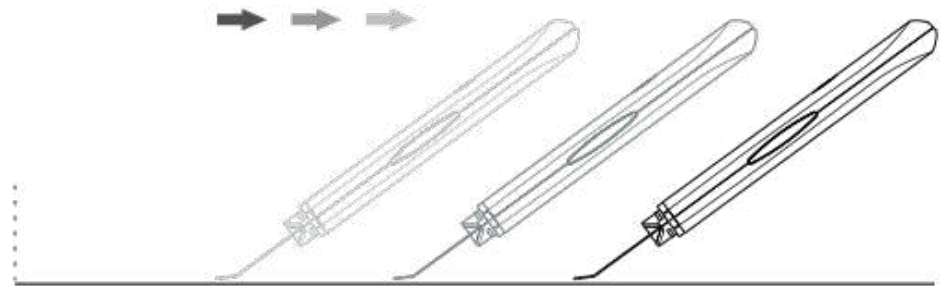
The procedure is performed using Medica+ Conductive Lifting Gel, serums, creams, essences, and facial masks.

In this mode, the scrubber can be used 3-4 times a week.

## Mode 3 - Lifting Charge - Lifting Mode

This mode is used for skin lifting and tightening, wrinkle smoothing, and deeper penetration of skincare products into the skin layers.

In this mode, hold the scrubber with the flat side of the blade against the skin. Move the scrubber slowly from the center to the periphery.

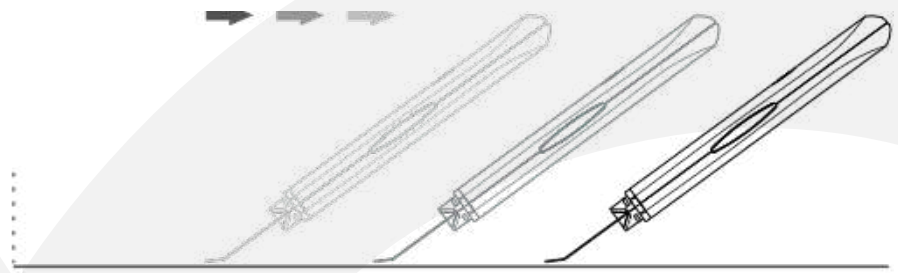The procedure is performed using serums, creams, essences, and face masks.

In this mode, the scrubber can be used daily.

# Technical Specifications

Weight – 0.22 kg

Output Voltage – 5V, 1A

Device Dimensions (H x W x D) – 18 x 5.3 x 1.6 cm

Package Dimensions (H x W x D) – 20.5 x 9.5 x 5.2 cm

Frequency – 23 kHz

Power – 3W

Battery – 500 mAh

# Warranty Terms

Dear customers!

Please keep this manual properly to use the warranty card. The warranty is provided for one year from the date of purchase (or two or three years if an extended warranty is purchased) in case of malfunction under normal operating conditions.

The warranty does not cover the following cases:

1. If the product is damaged due to improper use or lack of proper maintenance.
2. If the device is damaged as a result of disassembly or modifications.
3. If non-original parts have been used with the device.
4. Damage caused by force majeure events.

According to the terms, technical service provided outside the warranty period is subject to additional charges.

Any parts replaced during the warranty period remain the property of MEDICA+.

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