ENG / UA / PL

**Ultrasonic
Toothbrush**

**Probrush 10.0XN
Portable**

MEDICA+ LUX

JAPAN TECHNOLOGY

MEDICA+ LUX

JAPAN TECHNOLOGY

ENG

MEDICA+ 

JAPAN TECHNOLOGY

**Instruction Ultrasonic
Toothbrush**

Probrush 10.0XN Portable

Content

Introduction.....	01
Important tips.....	02-04
Reminder, Other recommends.....	05
Other functions indicate	06-07
List of components	08
How to use	09-10
Battery status and charging.....	10-11
Product specifications.....	12
Fault diagnosis	13-14
Warranty terms.....	15-16

Introduction

- Thank you for choosing **Probrush 10.0XN Portable** electric toothbrush! Please read the safety instructions before using this product. The instructions will help you understand the toothbrush functions and related warnings.

Important tips

- For users of orthodontics, tooth repair, dental implants, etc. Please use the toothbrush under guidance of the dentist.
- Any questions involving pacemakers or other implanted devices, con-sult your doctor or the manufacturer of the implant before use.
- This product is designed for oral cleaning only. Do not use for any other purpose. If you feel uncomfortable or pain, please stop using this product and contact your doctor.
- This product is not suitable for mental patients and children, unless it is used under the supervision of a guardian.

- Do not immerse the charging base in water or use it in any wet environment. Ensure that the charging base and the power plug are completely dry when access to power. If there is any damage, do not use .
- This product contains a build-in rechargeable battery, please do not get close to open fire and do not use or place it at a high temperature.
- If the product is damaged due to collision or fall, please stop using the product.
- Do not modify, disassemble or repair the product on your own.
- Please check if there is any obvious damage before using this product. If there is damage or serious fall, do not use, and contact the manufacturer back to repair.

- Before connecting the power supply, please check whether the local power supply voltage is consistent with the voltage marks in this product.
- Do not use this product while bathing. The toothbrush is waterproof, however, the steam produced may penetrate into the toothbrush and cause a short-circuit.
- Do not wash the charger base, do not reach for the charger that has fallen into the water, doing so may result in electrocution.

Reminder

- This product is used for personal cleaning only. It is not recommended for more than one person.
- This product should not be exposed to sunlight and should be placed in a dry and ventilated environment.

Other recommendations

- Professional dentist recommends brushing for 2-3 minutes at a time
- Professional dentist recommends sonic brushing method to clean the teeth.
- Professional dentist recommends replacing the brushhead every 2-3 months.

Other functions indicate

- It divides the oral cavity into four areas and jitters them every 30 seconds to remind you to change the area
- After reaching the recommended minimum brushing time, it will automatically shut down after two minutes

Pink light indicates: Basic oral cleaning, suitable for gingival sensitive people and initial users


Blue light indicates: Deep oral cleaning, suitable for most people

Green light indicates: Deep cleaning needs, two frequency vibration alternately, effectively remove dental plaque

Battery indicator: Light shows long red light when the battery is low; During charging the indicator light shows red light flashing state; When fully charged, the green light is always on

Charging instruction: Fully charged in about two hours

List of components

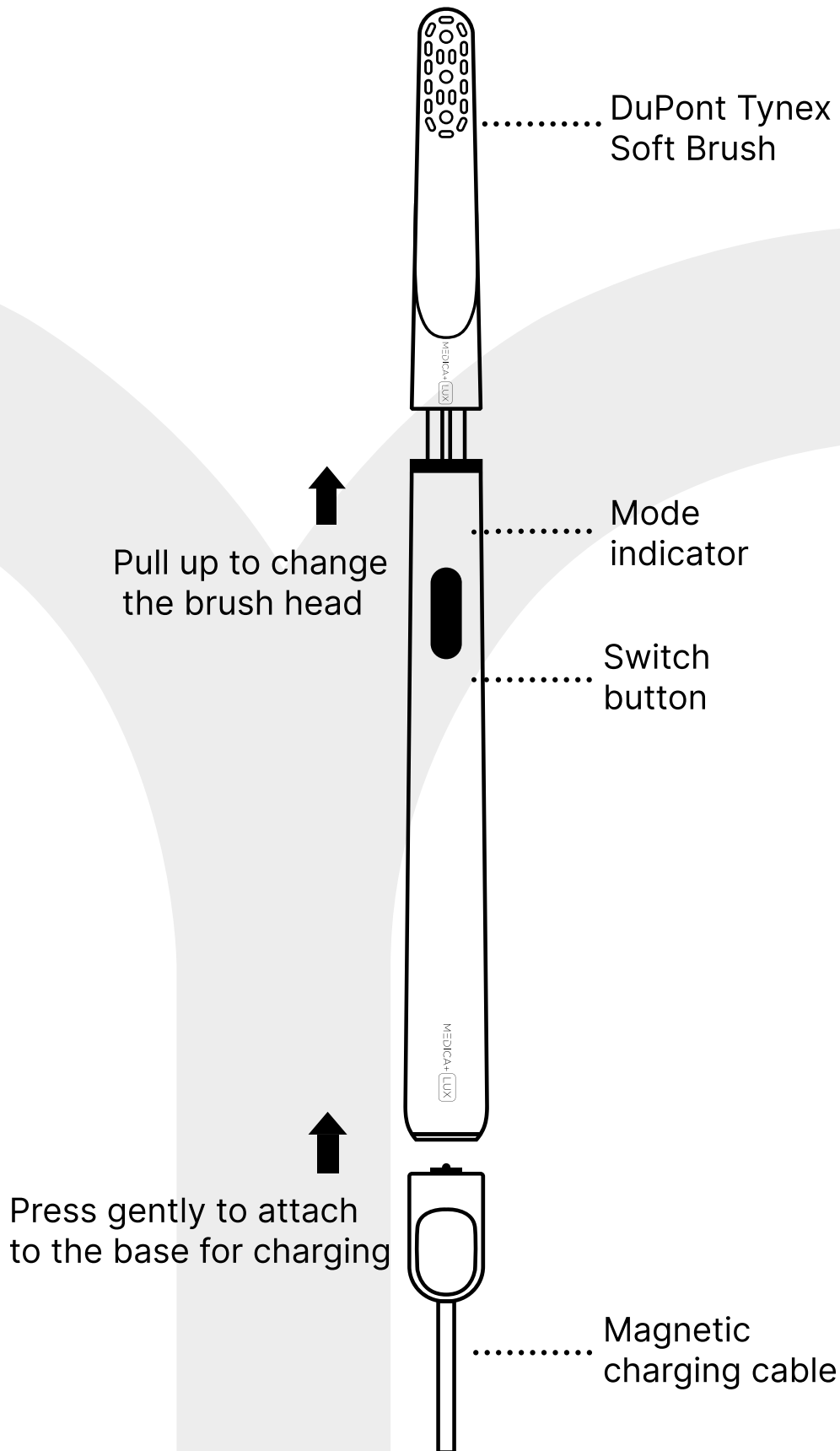How to use

1. Soak the bristles and squeeze a proper amount of toothpaste on the toothbrush.
2. Press the powerbutton to start the brush handle and select the desired mode.
3. Place the bristles close to the teeth, move the toothbrush slowly, and brush the teeth. Every 30 seconds during the process, the toothbrush shakes to prompt to change the brushing area, and it will automatically shut down after 2 minutes. If you need to continue brushing, press the powerbutton again.
4. To prevent toothpaste foam from spattering, turn off the toothbrush before removing from mouth.

5. After each use, rinse the toothbrush with clean water and wipe off with dry cloth.

*The bottom of this product is a nickel-metal hydride. Long term contact with salted toothpaste residue may cause the accessory to rust. Use a dry cloth to wipe off the end off metal part after use to prevent rust.

Battery status and charging

- This product has undergone a complete charge and discharge aging before leaving the factory, but due to the production cycle and natural discharge, the first time the battery endurance time used may be different.

- This product uses a USB interface for charging. You can use a charging adapter that meets national standards to charge it, but using a lower standard charger may bring unpredictable safety risks.
- According to the test, it takes about 2-3 hours to finish a complete charge. The specific charging time may be affected by the remaining power and the difference of charging environment.
- When the battery indicator flashes in red, it means that the battery power reaches close to the limit. Please charge the toothbrush in time to prevent the built-in lithium battery from malfunctioning and damage to extensive discharge.
- In a laboratory environment and when the battery is fully charged, based on the average use time of the three gears, 2 times a day, 2 minutes each time, the toothbrush battery can be used for 45-60 days.

Product specifications

Product name	Probrush 10.0XN Portable
Material	Aluminum alloy, ABS
Bristle	Dupont Tynex
Power	1W
Input	5V-0.5A
Battery capacity	320 mAh
Frequency	18000-22000 times/minute
Charging	USB Magnetic induction charging
Mode	Three Gears Modes
Water Proof	IPX 7
Weight	Toothbrush holder: 5 grms Toothbrush: 51 grms, Charger: 15.3 grms

Fault diagnosis

Problem	Possible cause	Remedy
Vibration intensity weakened	Battery is low	Charge the toothbrush for around 2 hours
After charging, the vibration is still weak	Falling or bumping causes the motor shifted	Send the product back for repair (do not disassemble it by your own so as not to cause secondary damage)
	The product is not placed vertically on the charger and fails to charge	Put the product upright so that its bottom is in contact with the charger base, and confirm that the indicator light is on
	Toothpaste or water stains on the bottom of the product or on the charger	Wipe clean the bottom and bottom contact surface of the product and try again

Fault diagnosis

Problem	Possible cause	Remedy
Toothbrush is not working	Just purchased or not used for a long time	Charge the product for at least 2 hours
Abnormal brush head noise	Brush head is incorrectly assembled	After removing it, put it back in place properly

Warranty Terms

Dear customers,

Please keep this instruction manual properly to use the warranty card. The warranty is provided for one year from the date of purchase (or for 2 or 3 years if the extended warranty is purchased) in case of malfunction under normal operating conditions.

The warranty does not cover the following cases:

1. If the product is damaged due to improper use or lack of proper maintenance.
2. If the product is damaged as a result of disassembly or modifications.

3. If parts not supplied by our company were applied to and used with the device.
4. Damage caused by force majeure circumstances.

According to the terms, technical maintenance provided outside the warranty is subject to additional charges.

Parts of the product replaced during the warranty period remain the property of MEDICA+.

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