
ENG / UA / PL

LEDMask 8X

medica+
JAPAN TECHNOLOGY

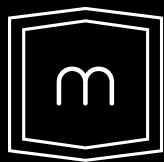medica+

JAPAN TECHNOLOGY

ENG

medica+
JAPAN TECHNOLOGY

Instruction manual

LED mask for light therapy

Content

Operating Principle	01-02
Equipment Contains	02
Light Spectrums	03-04
Usage Instructions.....	05-06
Contraindications	06
Precautions	07
Warranty Conditions	08
Technical Specifications	09

Dear Customer,
Thank you for choosing Medica+ products. Please read this manual carefully before using the device.

Operating Principle

The LED mask uses light therapy — a method that applies different wavelengths of LED light to the skin. Each light color penetrates the skin at a specific depth and activates natural regeneration processes. It is suitable for all skin types.


When the skin absorbs the light, it is converted into energy that “awakens” the cells:

- Stimulates the production of collagen and elastin (controlled by fibroblasts)
- Improves blood circulation and cellular nourishment
- Reduces pigmentation and evens out skin tone
- Leaves the skin firmer, smoother, and rejuvenated
- Promotes even melanin distribution

LED therapy is painless, does not cause heat or irritation, and requires no recovery time. It does not interfere with the body's natural functions — it simply activates existing natural processes in the skin.

- Initial results are usually visible after just a few sessions, and with regular use, more noticeable improvements typically appear within 3–4 weeks.

Equipment Contains

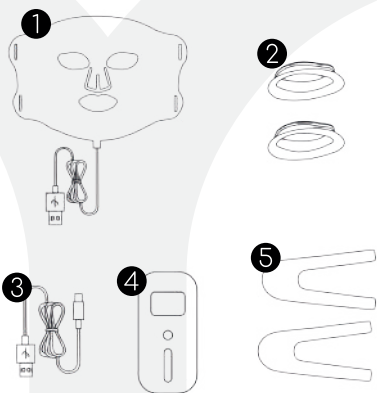1. LED light therapy mask with connecting cable
2. Eye inserts
3. Charging cable
4. Control module
5. Adjustable straps

Light Spectrums

● Red Light (630 nm)

Main effect: Stimulates fibroblasts and activates collagen and elastin production

Best for: Skin rejuvenation, reducing wrinkles, improving skin firmness and tone

● Green Light (520 nm)

Main effect: Evens out skin tone, brightens pigmentation, soothes the skin

Best for: Pigmented skin, post-sun exposure, anti-aging care

● Blue Light (460 nm)

Main effect: Antibacterial effect, reduces inflammation

Best for: Acne, breakouts, oily and problematic skin

Turquoise Light (500 nm)

Main effect: Hydrates and soothes the skin, reduces irritation

Best for: Dry, sensitive, dehydrated skin; post-treatment care

Yellow Light (590 nm)

Main effect: Improves lymphatic flow, reduces puffiness, stimulates detoxification

Best for: Tired skin, puffiness, oily and problematic skin

Purple Light (425 nm)

Main effect: A combination of red and blue light; promotes regeneration

Best for: Mature, weakened, or problematic skin; post-acne recovery

White Light (470 nm)

Main effect: Calms, relaxes, and enhances cellular activity

Best for: Stressed, sensitive skin; enhances the effect of skincare products

Usage Instructions

1. Cleanse your facial skin and apply toner.
2. We recommend applying a serum, fluid, or active care product according to your skin type — this will enhance the effect of the treatment.
3. Take out the mask and connect the USB cable to the controller.
4. Turn on the mask by pressing the P/B button.
5. Press the color selection button to choose the light mode, then set the session time (recommended: 15–20 minutes). The procedure should be performed 3–4 times a week.

6. Do not look directly at the light while the mask is operating.
7. A short press on the P/B button allows you to adjust the light intensity.
8. Perform the procedure while lying on a bed or sofa for complete relaxation.
9. After the session, remove the mask and wipe it with a damp antiseptic wipe or a cotton pad moistened with chlorhexidine. Do not rinse the mask under running water.

Contraindications

- Pregnancy
- Epilepsy and seizure disorders
- Acute inflammatory conditions or open wounds on the skin
- Severe skin diseases (eczema, psoriasis in the acute phase)
- Individual intolerance or allergy to the mask materials
- Use of photosensitizing medications
- Chronic diseases affecting the skin (doctor consultation required)
- Oncological diseases

Precautions

- Keep your eyes closed during the session. If the light feels too bright, cover your eyes with cotton pads.
- A slight silicone smell may be noticeable after opening the package — this is normal for a new mask and will quickly dissipate.
- If you use the mask in bed before sleeping, be sure to remove it and place it back in the box after the session to avoid damage.
- Do not soak the mask in water. To clean it, simply wipe it with antiseptic wipes or a cotton pad soaked in chlorhexidine.
- Do not use harsh detergents or chemical agents to clean the device.
- Avoid prolonged exposure of the mask to direct sunlight.

Warranty Conditions

Dear customers!

Please keep this manual carefully to use your warranty card.

The product is covered by a one-year warranty from the date of purchase (or for 2 or 3 years if an extended warranty was purchased), provided the device was used under normal operating conditions.

The warranty does not cover the following cases:

1. Damage caused by improper use or lack of proper maintenance.
2. Damage resulting from disassembly or modifications.
3. Use of parts or accessories not provided by our company.
4. Damage caused by force majeure events.

According to the terms, any service performed outside the warranty period will be subject to additional charges.

Parts replaced under warranty remain the property of MEDICA+.

Technical Specifications

- Product Name: LED Light Therapy Mask
MEDICA+ LEDMask 8X
- Number of LEDs: 108
- Material: Silicone
- Rated Input: DC 5V/1A
- Rated Power: 5W
- Model: E100
- Connector Type: Type-C
- Battery Capacity: 2000 mAh
- Product Dimensions: 295×218×5.3 mm

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