ENG / UA / PL

Portable Oral irrigator


Water Puls 8.0

medica+
JAPAN TECHNOLOGY

ENG

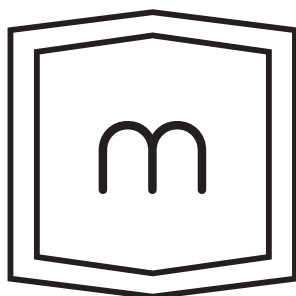medicat+
JAPAN TECHNOLOGY

Instruction

Water Puls 8.0

Content

Safety Warnings	01-02
Product Parameters	02
Component Name	03-04
Use Guide	05-09
Daily Maintenance	09
Disposal	10
Attention	10
Product Configuration	11


Safety Warnings

1. This appliance cannot be used by children and people with reduced physical, sensory or mental capabilities or lack of experience and knowledge unless they have been given or instructed concerning use of the appliance in a safe way and understand the hazards involved.
2. Children should be supervised to ensure that they do not play with the appliance. Cleaning and maintenance shall not be performed by children without supervision.
3. Do not use the charger outdoors or near heated surfaces.
4. This product is suitable for healthy people over 16 years old.
5. Do not clean the product or charger in the dishwasher.
6. Do not use the product for any purposes other than oral cleaning. Otherwise, accidents or injuries may be caused when water flows to the eyes, nose, ears or throat. Stop using the appliance and consult your doctor for support if you experience and pain or discomfort.
7. Do not share the nozzle with family or others. Otherwise, it may cause infection or inflammation.

8. Do not fill the water tank with water with a temperature of more than 40°C. Otherwise, it may cause burns.
9. Do not exceed the rating of the domestic power socket or wiring when charging
10. If you have had oral or gum surgery in the previous 2 months, consult your dentist before using this device.
11. For this product special protective measures for electromagnetic compatibility (EMC) shall be taken. If you have a heart rate regulator or other device, please consult your physician or the manufacturer of the device before use.

Product Parameters

Product Size	86*72*218mm
Input Voltage	5V \equiv
Input Power	2.5W
Product Material	ABS, PC, PMMA
Charging Time	\leq 8h
Waterproof Level	IPX7
Working Temperature	5-40°C

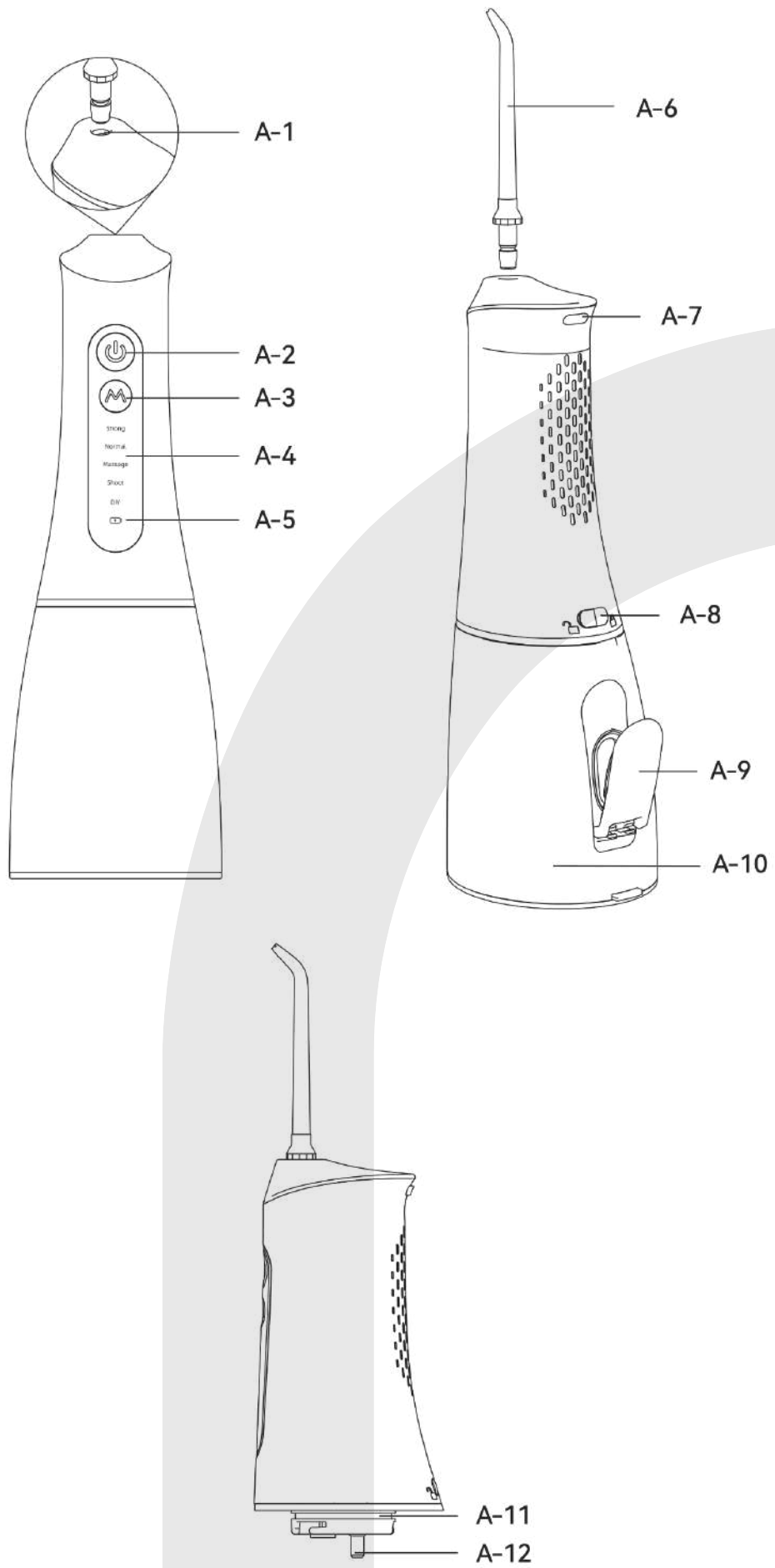Component Name

A-1 - Tip Hole

A-2 - On / off

A-3 - Mode Button

A-4 - Mode Indicator

A-5 - Charging Indicator

A-6 - The Tip

A-7 - Tips Release Button

A-8 - Outlet

A-9 - The Tank Cover

A-10 - The Water Tank

A-11 - Suction Pipe

A-12 - Filter

Use Guide

Before use please follow the steps below to charge the product:

Before use please follow the steps below to charge the product:





Product application part is BF type. Check and make the product and hands are dry, then connect the USB charging cable to the power supply adapter and plug the power plug into the power supply socket. Plug in the charging plug of the USB charging cable, and ensure plug in place.

TIPS

Red indicator light is found after use or the product can't be turned on properly. At this point, the battery has been insufficient, need to charge the product.

- It takes 7-8 hours to fully charge.
- When charging, the red charging indicator light is on.
- The green charging indicator light is on, indicating that the charging is complete
- When the battery is fully charged, the charging will be turned off for battery overcharge protection.

When using, please follow the steps:

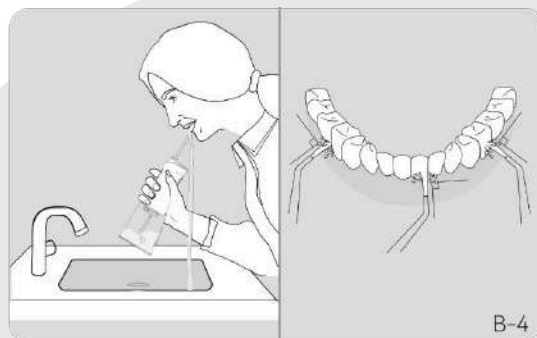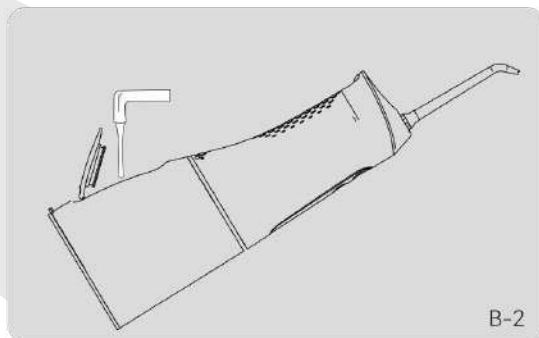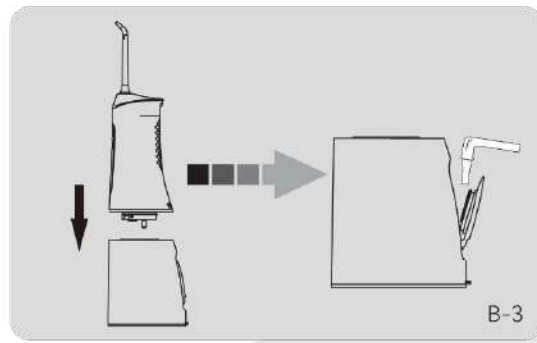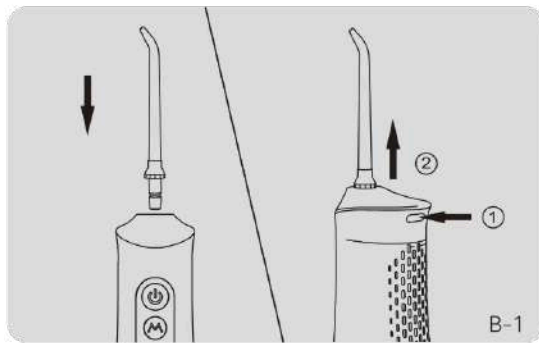1. Insert the tips into the hole at the top of the handle. If the tips need to be replaced, shut down the tips first and then press the release button of the tips lightly to easily replace the tips (B-1).
2. Put water into the tank. There are two ways to fill water, as follows: the first method: open the tank cover hold the tooth punch horizontally, put water into the tank and then cover the tank (B-2). Second method: remove the water tank hold the tank vertically, and load the water into the tank (B-3).

3. Includes five modes: Strong/Normal/Massage/ Shoot/DIY
- First mode(Strong): spray water three times per second, first indicator light, timing 2 minutes. Strong flow is for deep cleaning
 - Second mode(Normal): continuous water spraying, second indicator light, timing 2 minutes. Easily remove plaque for teeth, this setting can fulfill your daily life use
 - Third mode(Massage): spray water three times per second, third indicator light, timing 2 minutes. This mode is for massaging gums
 - Fourth mode(Shoot): when you choose this mode, fourth indicator light, you should press the On/off Button all the time, then it will be working. It focuses on cleaning stubborn stains.
 - Fifth mode(DIY): customize the water pressure from 80psi-150psi as your needs, fifth indicator light, timing 2 minutes. Choose the most suitable water pressure for yourself.

TIPS

If you use our product for the first time please try soft mode for a week until adaptation. Then try other modes.

4. Before use, hold the tooth punch vertically, with the tips sticking into the mouth and aiming at the teeth or gums. At the same time, open the mouth slightly so that water can flow out smoothly (B-4).
5. Press the switch button to use. Please note the following points when using:
 - align the gum line with the tips and be perpendicular to the gum,
 - move slowly along the teeth and clean them in sequence;
 - clean teeth, trays, crowns and dentures thoroughly.

TIPS

Do not spray water directly down the periodontal pouch. After use, please press the switch button to turn off the machine. This product has 2 minute timing function, use 2 minute automatic shutdown. If you still need to use it, please press the switch button. The product has the mode memory function.

After use, please follow the following steps:

1. Open the tank cap and pour the water out.
2. Press the switch button, the tooth punch works, and automatically drain the water inside.
3. Press the switch button to turn off the tooth punch and then close the tank cover.
4. Dry the tooth punch with a cloth.

TIPS

To avoid the spread of bacteria, make sure there is no residual water in the tooth punch. If you plan not to use the product for a long time, please stock the product after dry.

Daily Maintenance

1. Please clean the product before using. Or cleaning the product every 6 month if you don't use it for a long time.
2. Clean the product with clean water or neutral detergent. Do not use corrosive or abrasive detergent. Do not use hot water over 40°C in order to avoid damaging the product.

* Replace the tips every 6 months.

Disposal



1. The power source of this product is a Lithium-ion battery which can be recycled. At the end of the products life, please remove the battery before you discard the product, recycle or dispose the battery according to the local waste management requirements. If you have recycling questions, contact your local waste management office to ask help.
2. Before you decide to remove the battery, you need to prepare two tools: pliers and scissors.
3. When using the tools, please observe the basic safety rules and take protective measures to avoid injury.

Attention

Properly recycle or dispose of the battery. Ask your local waste management department whether to recycle other parts of the product.


Be careful not to short the positive and negative terminals of the battery. Don't charge the product after removing the battery.

Product Configuration

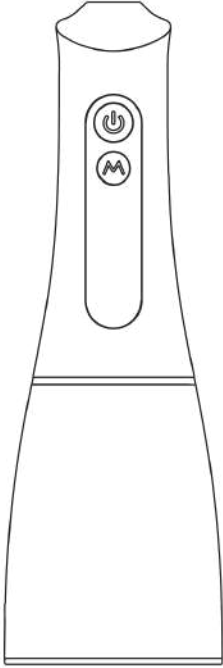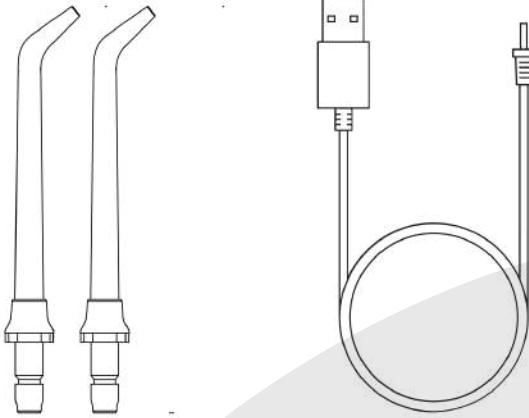	
Water flosser × 1 Charging cable × 1	Tips × 2

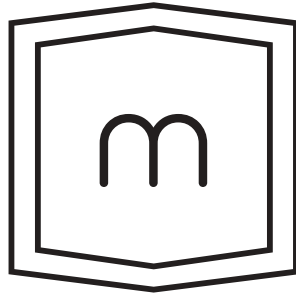1. Standard tip: Suitable for oral rinse food residue

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