
ENG / UA / PL


Portable darsonval

DARSOLINE 7.0

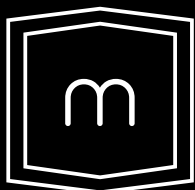medica+
JAPAN TECHNOLOGY

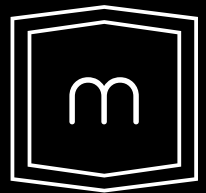medica+

JAPAN TECHNOLOGY

ENG

medicat+

JAPAN TECHNOLOGY

Instructions

DARSOLINE 7.0

Content

Purpose	01
Operating Principles.....	02
Application.....	03
General Recommendations	04-05
Operating procedure	06
Indications	07
Contraindications	08
How to perform the procedure.....	08-09
Important safety measures	10
Parameters	11
Package contents	12
Electrode "Spiral".....	13
Electrode "Roller"	14
Electrode "Arc"	15
Electrode "Petal"	16
Electrode "Comb".....	17
Electrode "Mushroom"	18
Electrode "Drop"	19
Warranty Terms.....	20

Dear Customer,
Thank you for choosing Medica+ products.
Please read this instruction carefully before using
the device.

PURPOSE

The device affects the skin with a high-frequency discharge without damaging the biological structures of the tissues.

It is designed for use in cosmetology, medicine, and home care for the prevention and treatment of dermatological, therapeutic, neurological, surgical, dental, and ENT diseases, as well as for hair health improvement and skin rejuvenation.

OPERATING PRINCIPLES

The device generates sequential electro-ionizing discharges with a pulse frequency of 50-60 kHz. The output voltage amplitude is smoothly adjusted using a regulator.

The therapeutic effect of the device on biological objects is based on the discharge that occurs between the body's surface and the electrode.

The device has no galvanic connection between the 220V power supply and the electrode output, meaning electric shock through the electrode is impossible.

APPLICATION

The method of using the device is divided into contact and non-contact approaches, with execution being either static or dynamic.

Contact Method

Place the electrode directly against the skin and move it smoothly without pressure across the surface, or hold it in place for targeted cauterization.

Non-Contact Method

Position the electrode 5-10 mm above the skin's surface and apply treatment using either a moving or static approach. The distance is determined by the discharge power set with the regulator.

Static Method

Involves targeted cauterization of the affected skin area by holding the electrode in place for a specific period.

Dynamic Method

Move the electrode using linear, circular, or "snake-like" motions.

General Recommendations for Using the Device

- Adjust the discharge power according to your individual perception, ensuring a mild tingling sensation on the skin that does not cause discomfort. Gradually, the power can be increased to a tolerable tingling level.
- Direct electrode contact with the mucous membrane of the eyes is strictly prohibited. However, application to the skin of closed eyelids is allowed at a low discharge power.
- The average duration of the procedure is 5-10 minutes. A standard treatment course consists of 10-15 sessions, conducted daily or every other day. If necessary, repeat the treatment course after one month.

Recommendation:

Gradually adjust the duration of exposure, movement speed, and discharge power based on personal sensations, modifying these parameters as needed. The most effective Darsonvalization procedures are those performed at the edge of comfort but without pain.

If redness appears on the skin, reduce the power or duration or increase the movement speed. Do not keep the electrode in one spot for too long, as this may cause redness and a burning sensation, except in the static procedure for targeted cauterization of the affected area.

PROCEDURE FOR WORKING WITH THE DEVICE

1. Insert the selected electrode fully into the socket located at the top of the device.
2. Plug the device into a power outlet.
3. Turn the regulator clockwise to switch on the device.
4. Adjust the required power level based on personal sensations of mild tingling. The device is now ready for use.

To turn off the device, turn the regulator counterclockwise until it clicks, then unplug it from the power outlet.

After turning off, you can remove the electrode from the socket or replace it with another one.

INDICATIONS

- Presence of acne
- Skin and subcutaneous inflammations (comedones)
- Enlarged pores
- Oily and problematic skin
- Pale and dull facial skin
- Post-acne marks
- Presence of scars and stretch marks on the body
- Loss of skin tone and elasticity
- Fine lines and age-related wrinkles
- Varicose veins
- Cellulite and stretch marks
- Dandruff, seborrhea on the scalp
- Hair loss
- Hair brittleness
- Headaches
- Insomnia
- Presence of bruises and hematomas on the body

CONTRAINDICATIONS

- Pregnancy
- Children under 6 years old
- Presence of oncological diseases
- Hypertension
- Blood clotting disorders or thrombophlebitis
- Cardiovascular diseases
- Presence of a pacemaker
- Epilepsy
- Tuberculosis
- Individual contraindications to the procedure (it is important to consult a doctor before the procedure)

HOW TO PERFORM THE PROCEDURE

1. Before the procedure, clean and dry the area where the darsonvalization will be performed.
2. Remove all metal objects (hairpins, chains, earrings, etc.).
3. Select the necessary attachment.
4. When using Darsonval on the face or body, it is recommended to apply Medica+ Antibacterial Suspension beforehand, which has wound-healing and regenerative properties and contributes to a more comfortable and pleasant darsonvalization procedure.

5. When using it on the scalp, brush your hair.
6. When using it on the thigh, arms, or legs, no additional products are needed.
7. Start the procedure on the problem area for 5-10 minutes (increase the duration gradually). The duration of exposure in the neck area should not exceed 5 minutes. All device movements should be performed only along massage lines and on the problem area.
8. After completing the procedure on the face, it is recommended to apply Medica+ Soothing Mask or Medica+ Pore Control Mask.
9. After completing the procedure on the scalp, you can apply a hair-strengthening product.
10. After completing the procedure on the thigh area, you can apply an anti-cellulite or nourishing cream.

IMPORTANT SAFETY MEASURES

1. Do not use the device if you notice damage to the glass electrode.
2. Before storing the device, make sure it is turned off. Store the device in a dry place out of reach of children.
3. Use antiseptic wipes or a cotton pad soaked in chlorhexidine solution to clean the glass electrodes. Do not use corrosive agents for cleaning. Clean the device after each use.
4. Do not throw the device into fire or store it in high-temperature or humid environments.
5. Do not immerse the device in water or rinse it directly under a faucet.
6. If the device has been damaged due to a fall, do not use it.
7. Do not disassemble or attempt to repair this product yourself. For repairs, please contact the Medica+ service center.

PARAMETERS

Power supply voltage and frequency: 220V 50Hz

Input voltage: 0.71A

Power: 10W

Output signal form: Pulsed with high-frequency filling

Frequency: 50-60 Hz

Pulse filling frequency: 144 kHz

Pulse burst frequency: 100 Hz

Output voltage: 44 kV

Tip material: Korean ABS plastic

Electrode material: Glass

Package dimensions: 38×21×6 cm

Gross weight: 1 kg

Power supply: Mains-powered

PACKAGE CONTENTS:

1. Portable Darsonval Medica+ DARSOLINE 7.0 – 1 pc.
2. Spiral Electrode – 1 pc.
3. Roller Electrode – 1 pc.
4. Arc Electrode – 1 pc.
5. Petal Electrode – 1 pc.
6. Comb Electrode – 1 pc.
7. Mushroom Electrode – 1 pc.
8. Drop Electrode – 1 pc.

Electrode "Spiral"

- For the treatment of acne, pimples, and small rashes.
- For stimulating acupuncture points on the forehead.
- Quick treatment of the body and face.

Electrode "Roller"

- For the treatment of all body areas.
- Stimulates metabolism.
- Reduces swelling.

Electrode "Arc"

- Effective treatment of the neck, arms, and legs.
- Stimulates blood circulation.
- Improves well-being.

Electrode "Petal"

- Fights wrinkles.
- Evens out skin tone.
- Reduces redness.
- Suitable for treating narrow areas.

Electrode "Comb"

- Suitable for treating dandruff.
- Hair loss.
- Slow hair growth.

Electrode "Mushroom"

- Reduces inflammation.
- Removes acne and post-acne.
- Increases collagen production.
- Oxygenates the skin.

Electrode "Drop"

- Fights cellulite.
- Fights varicose veins.
- Relieves pain symptoms.
- For skin diseases.

WARRANTY TERMS:

Dear Customers!

Thank you for supporting the Medica+ brand.

Please keep this card properly to use the warranty service along with the warranty certificate.

The warranty is valid for one year from the date of purchase in case of malfunction under normal operating conditions.

The Warranty Does Not Cover the Following Cases:

1. If the product is damaged due to improper use or lack of proper maintenance.
2. If the product is damaged as a result of disassembly or modifications.
3. If non-original parts not supplied by our company have been used in the device.
4. Damage caused by force majeure circumstances.

According to the terms, any maintenance services provided outside the warranty are subject to additional charges.

All parts of the product replaced during the warranty period remain the property of Medica+.

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